

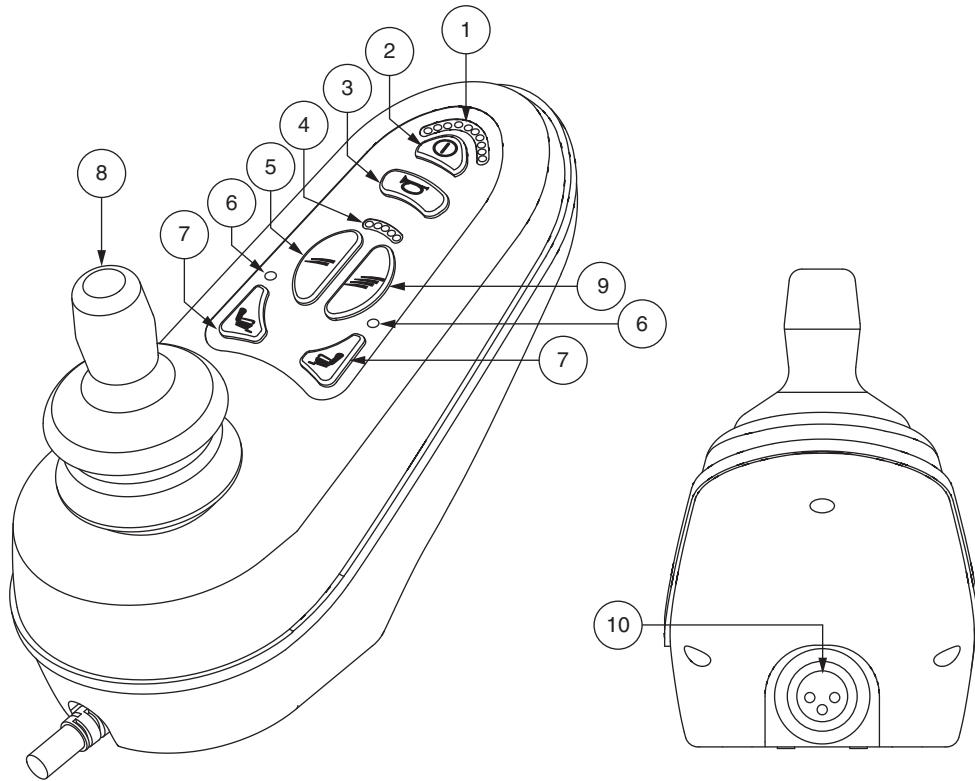
Product Insert


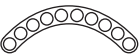









VR2



air[®]

Identification Key



 Controller	 Battery Condition Meter	 On/Off Button	 Horn Button
 Maximum Speed/ Profile Indicator	 Speed/Profile Decrease Button	 Mode LED	 Actuator Buttons
 Joystick	 Speed/Profile Increase Button	 Charger Port	



Product Insert

NOTE: *The following information is an addition to the information contained in the Basic Operating Instructions provided with your VR2 controller. It explains how to operate your Jazzy Air power elevating seat. Please retain this information for your records and contact your authorized Pride Provider with any questions.*

Jazzy Air

Your power chair is equipped with the feature that will be used in conjunction with your Jazzy Air. The Jazzy Air feature is mounted on a Pride power base for enhanced stability when the seating system is elevated.

There are certain situations, including some medical conditions, where the power chair user will need to practice operating the power chair equipped with Jazzy Air in the presence of a trained attendant. A trained attendant can be defined as a family member or care professional specially trained in assisting a power chair user in various daily living activities.

As you begin using Jazzy Air during daily activities, you will probably encounter situations in which you will need some practice. Simply take your time and you will soon be in full and confident control.

Elevation

The Jazzy Air feature can be activated by selecting the elevating seat actuator.

To select and adjust the seat elevation using the VR2 Controller:

1. Push the On/Off button once to power on the controller.
2. Push the Actuator Button until the Mode LEDs are illuminated.
3. When the elevating seat mode is selected, give a forward command to the joystick to adjust position in one direction or give a reverse command to the joystick to adjust position in the opposite direction.
4. If you want to return to another profile, push and release either Actuator Button until you return to the desired drive profile.

NOTE: *If your VR2 settings are different than those described, contact your authorized Pride Provider. Your authorized Pride Provider may have changed the programming of your VR2 Controller.*

To activate the Jazzy Air using a single switch:

1. Press the single switch once to raise the seat and activate Jazzy Air.
2. Press the single switch once to lower the seat.

NOTE: *The Jazzy Air feature is automatically engaged when the power elevating seat is activated, and an LED on the Maximum Speed/Profile Indicator is blinking on your controller.*

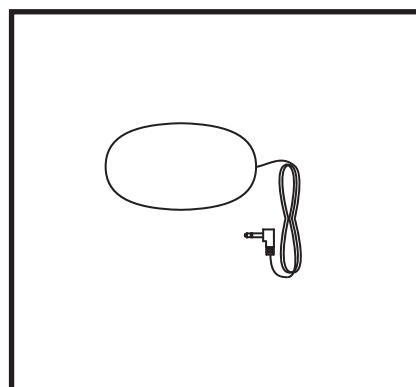


Figure 1. Single Switch

The power chair should be used on an indoor, level surface, such as malls, offices, and/or medical facilities, to activate the Jazzy Air. Once the seating system actuator is started, the Jazzy Air limits movement of the front caster wheels and restricts the power chair of climbing ability.

WARNING! The Jazzy Air feature should be used only on flat, level, indoor surface, such as in malls, offices, and/or medical facilities. Never raise the seat from its lowest position on an inclined surface. Failure to heed this warning can result in the power chair tipping over.

WARNING! The Jazzy Air feature is intended for limited outdoor use on only flat, level, hardscape surfaces. If you are unsure of your environment, do not use Jazzy Air outdoors.



PROHIBITED! Never use the Jazzy Air outdoors on uneven terrain (gravel, hills or other uneven, non-compact terrain) or in fields.



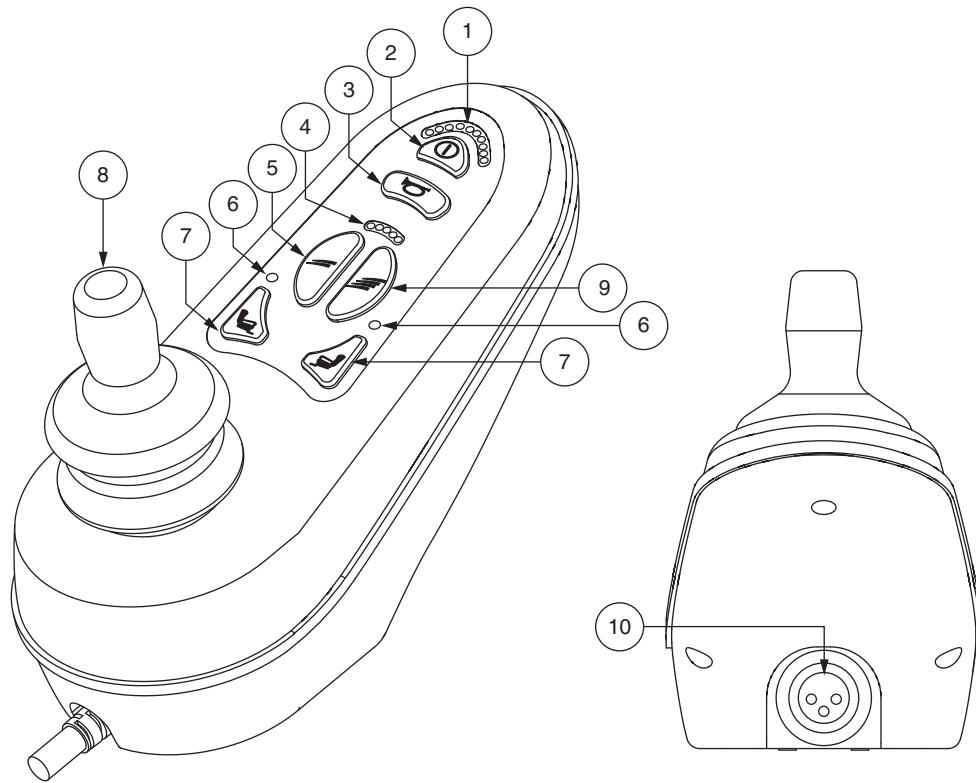
WARNING! Never raise the seat into Jazzy Air when operating your power chair on bumpy or uneven surfaces. Failure to heed this warning can result in the power chair tipping over.

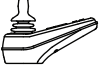










WARNING! Always fasten the positioning belt when operating the power chair.

WARNING! Do not operate the Jazzy Air or change the seat height in any way while your power chair is positioned under a fixed object, such as a table or a desk. Keep area clear prior to and during movement.

WARNING! When the Jazzy Air lockouts are not engaged, an audible alarm will notify you and the battery condition meter will scroll 10 bars. Speed will be reduced to 25%, when elevated. If this warning still appears on level ground, please discontinue use of Jazzy Air and consult your authorized Pride Provider for assistance.

NOTE: The Jazzy Air feature is equipped with a system that allows the power chair to drive at a predetermined, safe speed when the seat is elevated. Always check to be sure this system is operating properly before using your power chair.



 Controller	 Battery Condition Meter	 On/Off Button	 Horn Button
 Maximum Speed/ Profile Indicator	 Speed/Profile Decrease Button	 Mode LED	 Actuator Buttons
 Joystick	 Speed/Profile Increase Button	 Charger Port	

USA

401 York Avenue
Duryea, PA 18642

182 Susquehanna Avenue
Exeter, PA 18643-2694

Canada

5096 South Service Road
Beamsville, Ontario L0R 1B3

UK

32 Wedgwood Road
Bicester, Oxfordshire OX26 4UL
www.pride-mobility.co.uk

Australia

20-24 Apollo Drive
Hallam, Victoria 3803
www.pridemobility.com.au

New Zealand

38 Lansford Crescent
Avondale
Auckland, New Zealand 1007
www.pridemobility.co.nz

Italy

Via del Progresso, ang. Via del Lavoro
Loc. Prato della Corte
00065 Fiano Romano (RM)
www.pride-italia.it

B.V.

Castricummer Werf 26
1901 RW Castricum
The Netherlands
www.pride-mobility.nl

www.pridemobility.com



I N F I N F B 3 5 0 3