

Q6 edge®

Get the Edge.

With features like standard 4-pole motors, Mid-Wheel 6® drive design and ATX Suspension, the Q6 Edge® is engineered to meet the performance needs of the most active user. The Q6 Edge accepts a complete range of seating and electronics options to deliver exceptional rehab capability at an excellent value.



Q6 EDGE® shown with TRU-Balance® Power Positioning System and Q-Logic Drive Controls




THE POWER YOU TRUST

Specifications

DRIVE WHEELS	14" pneumatic or solid ⁵
LOW IMPACT OMNI-CASTERS	Front: 6" solid Rear: 6" solid
SUSPENSION	ATX Suspension (Active-Trac® with extra stability) with enhanced rear suspension
MAXIMUM SPEED ¹	Up to 5 mph / Up to 6 mph ⁵
BRAKES	Intelligent braking (electronic, regenerative disc brakes)
GROUND CLEARANCE	3" at center of frame
TURNING RADIUS ²	20"
OVERALL SIZE ²	Length: 35.5" Width: 24"
SEATING	Synergy® Seat Adjustable Width: 10"-22" Adjustable Depth: 10"-22" TRU-Balance® Power Positioning System Adjustable Width: 14"-22" Adjustable Depth: 14"-22" High-back Seat (22" x 22" max. dim.) Medium-back Seat (22" x 22" max. dim.)
DRIVETRAIN	2-motor, 4-pole, Mid-Wheel 6®
BATTERIES	Two 12 volt, deep cycle NF-22
RANGE PER CHARGE ^{1,3}	Up to 15.5 miles
BATTERY CHARGER	8A, off-board
AVAILABLE ELECTRONICS	75A Q-Logic NE 75A Q-Logic NE+ 75A Q-Logic Drive Control System
WEIGHT CAPACITY	300 lbs. maximum (user weight)
BASE WEIGHT	138 lbs.
BATTERY WEIGHT ⁴	38 lbs. (each)
WARRANTY	Lifetime limited warranty on frame; 2-year warranty on electronics; 18-month warranty on drive motor

Shroud Colors

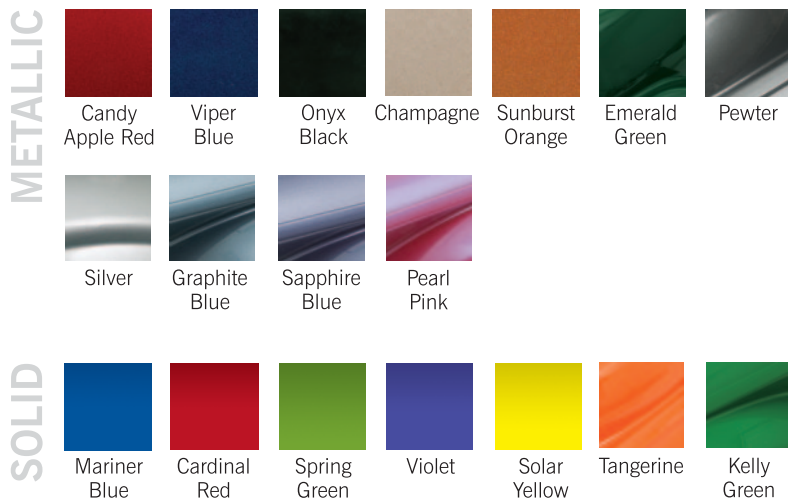

Features

- Mid-Wheel 6® allows six wheels on the ground for maximum stability
- Compatible with TRU-Balance® Power Positioning
- Low impact OMNI-Casters (nylon, spherical-shaped casters) on front and rear prevent wheel hang-ups
- ATX Suspension (Active-Trac® with extra stability) incorporates front OMNI-Casters and semi-independent rear caster beam for enhanced performance over more varied terrain

Options

- TRU-Balance® Lift & Tilt
- TRU-Comfort Seating System
- Q-Logic EX Drive Controls
- 6 mph motor package
- Accu-Trac Advanced Tracking Technology ⁶
- Power elevating seat
- Articulating vent tray ⁶
- Compact vent tray ⁶
- LED full lighting package
- 4" wide flat-free knobby extended traction tires ^{6,7}
- WC19 transit securement package ^{6,8}

Code Coverage

Q6 EDGE 3S - SS	(Code: K0848)
Q6 EDGE 3S - C	(Code: K0849)
Q6 EDGE 3SP - SS	(Code: K0856)
Q6 EDGE 3SP - C	(Code: K0857)
Q6 EDGE 3MP - SS	(Code: K0861)

FOOTNOTES:

- Range and speed vary with user weight, terrain type, battery charge, battery condition and tire pressure.
- Due to manufacturing tolerances and continual product improvement, this specification can be subject to variance of + or - 3 %. Dimensions listed are for power base only. Overall measurements will vary based on seating and accessory selections.
- Tested in accordance with ANSI/RESNA, WC Vol. 2, Section 4, standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Testing conducted at maximum weight capacities. This specification can be subject to a variance of +10%, -5%.
- Battery weight may vary based on manufacturer.
- Optional.
- Call for availability.
- Available as part of Community Use Package with 6 mph motors. Increases width by 1.5".
- Please refer to the Transit Securement Package supplement (part number INFMANU3658) for additional information.