

# Q1121

## Exceptional Maneuverability

The Quantum® 1121 delivers exceptional maneuverability, high-performance operation and superb rehab capability. Compatible with a wide range of rehab electronics and seating options like TRU-Balance® Power Positioning Systems, the 1121 is designed to meet the diverse needs of the active rehab client.




Quantum® 1121 shown with  
TRU-Balance® Power Positioning  
Systems and PG VSI controller

## Specifications

<b>DRIVE WHEELS</b>	14"
<b>REAR CASTERS</b>	8" solid
<b>ANTI-TIPS</b>	6" solid, adjustable, front-mounted
<b>SUSPENSION</b>	Active-Trac® Suspension
<b>MAXIMUM SPEED</b> <sup>1</sup>	Up to 5 mph
<b>BRAKES</b>	Intelligent Braking (electronic, regenerative disc brakes)
<b>GROUND CLEARANCE</b>	3.5"
<b>TURNING RADIUS</b> <sup>2</sup>	22.5"
<b>OVERALL SIZE</b> <sup>2</sup>	Length: 38.125" Width: 23.8"
<b>SEATING</b>	Synergy® Seat Adjustable Width: 10"-20" Adjustable Depth: 10"-20" TRU-Balance® Power Positioning System Adjustable Width: 14"-20" Adjustable Depth: 14"-20" Contoured seating (max. size: 20" x 22") Manual recliner Solid seat pan option
<b>DRIVETRAIN</b>	2 motor, mid-wheel drive
<b>BATTERIES</b>	Two 12 volt, NF-22
<b>RANGE PER CHARGE</b> <sup>1</sup>	Up to 15.5 miles
<b>BATTERY CHARGER</b>	8A, Off-board Charger
<b>AVAILABLE ELECTRONICS</b>	70A, PG VSI controller 70A, PG Remote Plus controller
<b>WEIGHT CAPACITY</b>	300 lbs. maximum (user weight)
<b>BASE WEIGHT</b>	140 lbs.
<b>BATTERY WEIGHT</b>	53.5 lbs. (each) <sup>3</sup>
<b>WARRANTY</b>	Lifetime limited warranty on frame 2-year warranty on electronics 18-month warranty on drive motor

## Shroud Colors

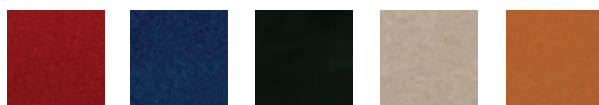

Candy Apple Red    Viper Blue    Onyx Black    Champagne    Sunburst Orange

(US) 866-800-2002 ▪ (Canada) 888-570-1113 ▪ [www.quantumrehab.com](http://www.quantumrehab.com)

The information contained herein is correct at the time of publication. We reserve the right to alter specifications without prior notice.

## Features

- 2-pole, patented two-motor, mid-wheel drive design provides superb power and acceleration
- Active-Trac® and rear caster suspension for outstanding outdoor performance
- Adjustable front anti-tips to meet your preferences and terrain challenges
- 8A, off-board charger
- Easy front access to batteries

## Options

- Cup holder
- Weather cover
- Swing-away joystick
- Elevating leg rests
- Oxygen tank holder
- Cane/crutch holder
- Walker holder
- Rear basket
- Angle-adjustable footplates
- Saddlebag

## Code Coverage

Q1121 3SP - SS	<b>(Code: K0856)</b>
Q1121 3SP - C	<b>(Code: K0857)</b>
Q1121 3MP - SS	<b>(Code: K0861)</b>

### FOOTNOTES:

- 1) Range and speed vary with user weight, terrain type, battery charge, battery condition and tire pressure.
- 2) Due to manufacturing tolerances and continual product improvement, this specification can be subject to variance of + or - 3 %. Dimensions listed are for power base only. Overall measurements will vary based on seating and accessory selections.
- 3) Battery weight may vary based on manufacturer.

**QUANTUM**<sup>®</sup>  
LIFE BEYOND LIMITS<sup>™</sup>