

# QUANTUM<sup>®</sup> DYNAMO<sup>®</sup>

## The Power Chair for Playing

Quantum<sup>®</sup> Rehab realizes the potential of the versatile Quantum product line to serve the needs of children with a variety of mobility issues. The new Quantum Dynamo<sup>®</sup> and Dynamo ATS are the first pediatric power bases to offer Pride's patented mid-wheel drive

design. Both the Dynamo and Dynamo ATS come standard with the pediatric Synergy<sup>™</sup> seating system, fully programmable electronics, pediatric swing-away footrests and the dependability of a Quantum product.



## Why choose the Quantum<sup>®</sup> Dynamo<sup>®</sup> for pediatrics?

The Quantum Dynamo is an ideal power base for pediatric use. Its compact design fits easily into the close quarters of a child's environment, while its 16.5" turning radius gives the unit smooth maneuverability

in the home, at school or on the playground. To keep children eye-to-eye with their friends, both the Dynamo and Dynamo ATS offer an optional seat-to-floor height as low as 15.5".

## The Quantum<sup>®</sup> Dynamo<sup>®</sup> ATS

The Quantum Dynamo ATS offers all of the same features as the Dynamo, with the added benefit of Active-Trac<sup>®</sup> Suspension, while maintaining an impressive 18" turning radius. Pride's Active-Trac<sup>®</sup> Suspension allows active children the opportunity to negotiate the roughest terrain to effortlessly keep up with their friends and family.



Quantum Dynamo ATS with Snug Seat Panda in Tilt/Recline Position



**Tiger Folding Mobility Frame**  
The frame folds for ease of storage and transport.

## Other features of the Dynamo and Dynamo ATS include:

- Stylish and sporty non-medical appearance.
- Seating systems as small as 10" x 10" that can grow with the child.
- Available with all of Quantum Rehab's specialty drive controls.
- Able to interface with a variety of power positioning systems including tilt, tilt/recline and tilt/recline with elevating seat combos.
- 10" knobby drive wheels and patented mid-wheel drive technology offer superior traction over most surfaces that today's active child wants to explore.
- Comes standard with a Synergy<sup>™</sup> seating system and swing-away footrests or drop-in footrests.

Quantum Dynamo ATS shown in T-Rex Red with Synergy<sup>™</sup> manual tilt seating system




Quantum Rehab is devoted solely to customizing the Quantum Power Chair line with specialty seating systems, controllers and a variety of rehab accessories based on an individual's needs and preferences. Our goal is to see that no child's needs go unmet.

For more information on the Quantum Dynamo and Dynamo ATS, call your **QUANTUM SPECIALIST** toll-free, 1-866-800-2002.

Quantum Dynamo ATS shown in Saber Yellow with Synergy<sup>™</sup> seating system and PG Remote Plus controller



- The Dynamo and Dynamo ATS come in a variety of fun colors.

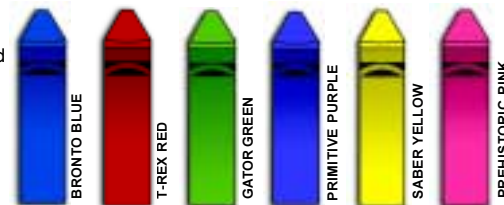

\*Colors may vary slightly from those shown.



# QUANTUM DYNAMO

THE POWER CHAIR FOR PLAYING

	DYNAMO®	DYNAMO® ATS
Wheels	10" pneumatic	10" pneumatic
Rear Casters	8" solid	8" solid
Anti-tips	6" solid, front-mounted	6" solid, front-mounted
Suspension	Limited	Active-Trac
Maximum Speed	Up to 4 mph	Up to 4 mph
Brakes	Intelligent Braking	Intelligent Braking
Ground Clearance	2"	2.5"
Turning Radius	16.5"	18"
Overall Size	Length: 32" Width: 23"	Length: 36.5" Width: 23"
Drive Train	Two-motor, mid-wheel drive	Two-motor, mid-wheel drive
Battery	Two 12-volt, 32AH(U-1)	Two 12-volt, 32AH(U-1)
Per-charge Range	Up to 25 miles	Up to 25 miles
Battery Charger	On-board, 4-amp	On-board, 4-amp
Standard Electronics	50-amp, PG VSI	70-amp, PG Remote Plus
Standard Seating	Synergy™ (14"W x 16"D)	Synergy™ (14"W x 16"D)
Seat to Floor Height	15.5"-17.5"	15.5"-17.5"
Weight Capacity	160 lbs.	160 lbs.
Transport Weight	Base: 79.5 lbs. Synergy™ Seat: 44 lbs.* (at 18"x18")	Base: 90.5 lbs. Synergy™ Seat: 44 lbs.* (at 18"x18")

\*Seat weight varies due to width, depth and height.  
The information contained herein is correct at the time of publication;  
we reserve the right to alter specifications without prior notice.

Speed & range vary with user weight, terrain type, battery charge, battery condition and tire pressure.

It is recommended that you consult a qualified Rehab Provider prior to fitting or ordering any Quantum Rehab product.

1-866-800-2002 (U.S.) • 1-888-570-1113 (Canada)  
[www.quantumrehab.com](http://www.quantumrehab.com)

DEALER STICKER

DYN-02-04



# QUANTUM DYNAMO

THE POWER CHAIR FOR PLAYING



**QUANTUM**  
**REHAB**  
Moving at the Speed of Life®

A Division of Pride Mobility Products® Corp.