


	Legend® 3-wheel	Legend® 4-wheel
Model number:	SC3000	SC3400
Weight capacity:	350 lbs.	350 lbs.
Maximum speed:	Up to 5.75 mph	Up to 5.75 mph
Ground clearance:	3.5"	3.5"
Turning radius:	41"	54"
Overall length:	45.5"	47"
Overall width:	24"	24"
Wheels Front: Rear: Anti-tips:	3.5" x 10" solid 4" x 10" solid Rear	3.5" x 10" solid 4" x 10" solid Rear
Suspension:	Rear	Rear
Transport weight Total w/o batteries:	157 lbs.	166 lbs.
Heaviest piece when disassembled:	61 lbs. (Rear section)	61 lbs. (Rear section)
Drive train:	Sealed transaxle, 24 volt DC motor	
Dual braking system:	Electronic, regenerative and electromechanical (electronic, regenerative disc brakes)	
Power elevating seat:	Option: Yes	
Std. Seating:	CRS deluxe reclining high-back with headrest and sliders	
Battery charger:	On-board, 3 amp (standard)	
Per-charge range:	Up to 22 miles	
Battery requirements:	(2) 12V deep cycle (AGM or gel cell type recommended)	
Battery size & weight: ¹	U-1; 25 lbs. each	
Warranty:	3-year limited warranty	

Colors:

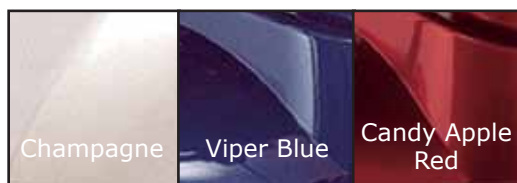

Benefits


- State-of-the-art rear suspension
- Top speed of up to 5.75 mph
- Maximum weight capacity of 350 lbs.
- Ultra-simple lockup for quick disassembly
- Advanced new reclining high-back seat with height-adjustable headrest
- Easily accessible tie down points (for transport of unoccupied scooter)
- Reduced length:
3-wheel: 45.5" / 4-wheel: 47"
- Dual-bulb, adjustable lower headlight for optimal pathway illumination
- Easy access to rear electronics
- Tiller-mounted, recessed "marker" headlight
- Tight turning radius:
3-wheel: 41" / 4-wheel: 54"
- Reinforced front nose section for durability

Options

- Power elevating seat
- Forearm crutch holder
- Cane/crutch holder
- Wishbone crutch holder
- Oxygen tank holder
- Rear basket
- Safety flag
- Walker holder
- Weather cover
- Saddlebag
- Lap belt

Silver Star Lifts & Ramps

- Backpacker® Plus
- Full Platform
- Silver Boom 300
- Backpacker® MV

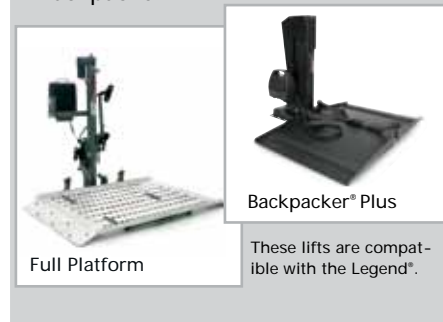

¹Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tire condition. This specification can be subject to a variance of (+ or -) 10%.

²Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or - 3%).

³AGM or Gel-Cell type required

⁴Battery weight may vary based on manufacturer.

⁵This range was obtained with the 4-wheel model.

⁶Tested in accordance with ANSI/RENSA, WC Vol 2, Section 4, and ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight capacity.

⁷Includes standard seat weight



The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed & range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3%.



EXPERIENCE THE LEGEND,
AND LIFE, LIKE NEVER BEFORE!

LEGEND®

The Legend® Scooter combines the latest style and performance technology at a great value. Standard features include a top speed of up to 5.75 mph, 350 lbs. maximum weight capacity and an easy lockup system for quick disassembly. Plus, the lower dual-bulb headlight offers excellent pathway illumination. The Legend's sleek styling and advanced features will have you turning heads every time you turn a corner.



Legend®

Pride®
Mobility Products Corp.

800-800-8586 (US) • 888-570-1113 (Canada) • www.pridemobility.com

LEGENDSS_08_09