

# QUANTUM<sup>®</sup>

1103 ULTRA



The versatile solution  
for everyday activities.

The Quantum<sup>®</sup> 1103 Ultra's compact, sporty mid-wheel drive design offers excellent tight quarter maneuvering and performance features like standard Active-Trac<sup>®</sup> Suspension provide enhanced outdoor capability. The Quantum 1103 Ultra delivers full rehab adaptability and stability on a very compact frame.


Quantum<sup>®</sup> 1103 Ultra shown with TRU-Balance<sup>®</sup>  
Power Positioning and PG Remote Plus controller

**QUANTUM<sup>®</sup>**  
LIFE BEYOND LIMITS<sup>™</sup>

## Specifications

|                                      |  |
|--------------------------------------|--|
| <b>DRIVE WHEELS</b>                  | 10"  |
| <b>REAR CASTERS</b>                  | 8" solid   |
| <b>ANTI-TIPS</b>                     | 5" solid, adjustable, front-mounted  |
| <b>SUSPENSION</b>                    | Active-Trac® Suspension  |
| <b>MAXIMUM SPEED</b> <sup>1</sup>    | Up to 4 mph (with high-speed motors)   |
| <b>BRAKES</b>                        | Intelligent braking<br>(electronic, regenerative disc brakes)  |
| <b>GROUND CLEARANCE</b> <sup>3</sup> | 2.5"   |
| <b>TURNING RADIUS</b> <sup>3</sup>   | 20.75"   |
| <b>OVERALL SIZE</b> <sup>3</sup>     | Length: 35.5" Width: 23"   |
| <b>SEATING</b>                       | Synergy® Seat<br>Adjustable Width: 10"-20"<br>Adjustable Depth: 10"-18"<br>TRU-Balance® Power Positioning System <sup>2</sup><br>Adjustable Width: 14"-20"<br>Adjustable Depth: 14"-20"<br>Contoured seating (20" x 20" max. dim.)<br>Manual recliner<br>Solid seat pan option |
| <b>DRIVETRAIN</b>                    | 2-motor, mid-wheel drive   |
| <b>BATTERIES</b>                     | Two 12 volt, deep cycle U-1  |
| <b>RANGE PER CHARGE</b> <sup>1</sup> | Up to 15 miles   |
| <b>BATTERY CHARGER</b>               | 5A, off-board  |
| <b>AVAILABLE ELECTRONICS</b>         | 50A, PG VSI<br>70A, PG Remote Plus   |
| <b>WEIGHT CAPACITY</b>               | 300 lbs. maximum (user weight)   |
| <b>BASE WEIGHT</b>                   | 92 lbs.  |
| <b>BATTERY WEIGHT</b> <sup>4</sup>   | 24.5 lbs. (each)   |
| <b>WARRANTY</b>                      | Lifetime limited warranty on frame;<br>2-year warranty on electronics;<br>18-month warranty on drive motor   |

## Shroud Colors

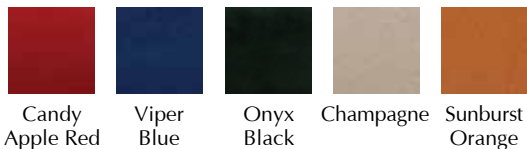

## Features

- 300 lbs. weight capacity
- Powerful two-motor, patented mid-wheel drive design offers outstanding maneuverability
- Low center of gravity design delivers a very smooth and stable ride
- 10" knobby drive wheels and optimum weight distribution provide superior traction
- Speed range adjustable up to 4 mph
- Standard lap belt
- Off-board battery charger
- Power elevating seat

## Options

- Cup holder
- Weather dust cover
- Swing-away joystick
- Swing-away footrests
- Elevating leg rests
- Oxygen tank holder
- Cane/crutch holder
- Walker holder
- Rear basket
- Angle-adjustable footplates
- Saddlebag
- Off-board charger (required with TRU-Balance® Power Positioning System)
- TRU-Balance® Power Positioning System available in conjunction with power elevating seat

### FOOTNOTES:

- 1) Range and speed vary with user weight, terrain type, battery charge, battery condition and tire pressure.
- 2) 200 lbs. weight capacity with TRU-Balance® Power Positioning
- 3) Due to manufacturing tolerances and continual product improvement, this specification can be subject to variance of + or - 3 %.
- 4) Battery weight may vary based on manufacturer