



Quantum® 6000

POWER CHAIR

Q6000 shown with TRU-Balance® 2 Power Tilt & Recline and Q-Logic drive controls

Available with high speed motors up to **8 mph**

Limitless rehab capability.


The aggressively styled, Quantum® 6000 incorporates Mid-Wheel® 6 design and a powerful drivetrain for excellent overall performance and maximum torque. The Quantum 6000 is fully compatible with a broad range of seating and electronic options for limitless rehab capability.



Specifications

DRIVE WHEELS	14" pneumatic/ XL: 14" solid
OMNI-CASTERS	Front: 5" solid Rear: 6" solid
SUSPENSION	Active-Trac® ATX Suspension with integrated rear dual-strut suspension and articulating caster beam
MAXIMUM SPEED	Up to 8 mph ^{1,2} / XL: Up to 6.5 mph ¹
BRAKES	Intelligent Braking (electronic, regenerative disc brakes)
GROUND CLEARANCE	3"
TURNING RADIUS³	22"
OVERALL SIZE³	Length: 37" Width: 25.875"
SEATING	Synergy® Seat Adjustable Width: 10"-22"/ XL: 10"-24" Adjustable Depth: 10"-20"/ XL: 10"-24" TRU-Balance® 2 Power Positioning System Adjustable Width: 14"-22"/XL: 14"-24" Adjustable Depth: 14"-20"/XL: 14"-24" High-back Seat (22" x 20")/ XL: (24" x 24")
DRIVE TRAIN	4-pole, Hammer Mid-Wheel 6®
BATTERIES	Two 12 volt, deep cycle Group 24
RANGE PER CHARGE¹	Up to 20 miles/ XL: Up to 15 miles
BATTERY CHARGER	8A, off-board
AVAILABLE ELECTRONICS	100A, Q-Logic
WEIGHT CAPACITY	300 lbs. maximum (user weight)/ XL: 400 lbs.
BASE WEIGHT	170 lbs.
BATTERY WEIGHT	53.5 lbs. (each)
WARRANTY	Lifetime limited warranty on frame 2-year warranty on electronics 18-month warranty on drive motor

Colors

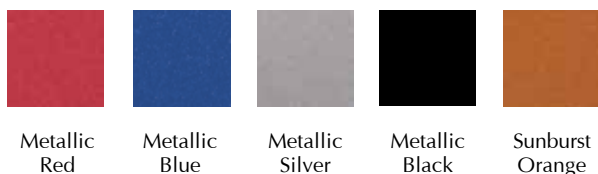

Features

- Mid-Wheel 6® allows six wheels on the ground for maximum stability
- Compatible with TRU-Balance® 2 Power Positioning
- OMNI-Casters (nylon, spherical-shaped casters) on front and rear prevent wheel hang-ups
- Side-mounted, easily accessible freewheel levers
- ATX Suspension (Active-Trac® with extra stability) incorporates front OMNI-Casters for enhanced performance over more varied terrain
- 4-pole Hammer motors
- Accepts (2) Group 24 & (1) NF-22 when used with a vent tray

Options

- Solid tires
- Swing-away joystick
- Multi-axis joystick mount
- Elevating leg rests
- Oxygen tank holder
- Neoprene leg rest sleeves
- Gel adductor buttons
- Vent tray
- TRU-Balance® 2 power elevating seat
- Angle-adjustable footplates
- High-mount foot platform
- 70° drop-in footrests
- TRU-Comfort™ seating system
- Contoured seating with solid seat pan
- Q-Logic Drive Control System

FOOTNOTES:

- 1) Range and speed vary with user weight, terrain type, battery charge, battery condition and tire pressure.
- 2) 8 mph with high speed motor package
- 3) Due to manufacturing tolerances and continual product improvement, this specification can be subject to variance of + or - 3 %.



(US) 866-800-2002
 (Canada) 888-570-1113
www.quantumrehab.com