




Hurricane® 4-wheel

Wrangler® 4-wheel

Model number:	PMV5001	PMV600
Weight capacity:	400 lbs.	400 lbs.
Maximum speed ¹ :	Up to 9 mph	Up to 10 mph
Ground clearance ² :	4"	7"
Turning radius:	65"	77"
Overall length ² :	52"	63.5"
Overall width ² :	24"	31"
Wheels Front: Rear: Anti-tips:	4" x 13" pneumatic 4" x 13" pneumatic Rear	6" x 13" pneumatic 6" x 13" pneumatic N/A
Suspension:	Front/Rear	Rear
Total weight w/o batteries ^{2,5} :	209 lbs.	216 lbs.
Heaviest piece when disassembled ² :	82 lbs. (front and rear section)	No frame disassembly
Drive train:	Sealed transaxle, 24 volt DC motor	2 in-line gear motors
Dual braking system:	Electronic, regenerative and electromechanical (electronic, regenerative disc brakes)	
Power elevating seat:	Yes (special order)	Yes (option)
Std. Seating:	CRS deluxe reclining high-back with headrest and sliders	Deluxe contoured high-back with headrest
Battery charger:	Off-board, 5 amp	On-board, 5 amp
Per-charge range ¹ :	Up to 15 miles with NF-22 Up to 10 miles with U-1	Up to 20 miles
Battery requirements ³ :	(2) 12V deep cycle (AGM or gel cell type recommended)	
Battery size & weight ⁴ :	NF-22; 38 lbs. each/ U-1: 24.5 lbs. each	Group 27; 63.5 lbs. each
Warranty:	3-year limited warranty	3-year limited warranty

Colors:

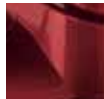

Candy Apple Red



Viper Blue

¹ Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition, and tire condition. This specification can be subject to a variance of (+ or -) 10%.

² Tested in accordance with ASNI/RESNA, MC Vol2, Section 4 & ISO 7176-4 standards. Results derived from theoretical calculation based on battery specifications and drive system performance. Test conducted at maximum weight capacity.

³ Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or -) 3%.

⁴ AGM or Gel-Cell type required. Battery weight may vary due to manufacturer.

⁵ Includes standard seat & weight.

Hurricane Benefits

- Delta tiller features wraparound handles for easy throttle operation with thumbs as well as fingers
- Four-pole motor and 110-amp controller for exceptional power, torque and efficiency
- Full lighting package including fender-mounted directional headlights
- Four-wheel design and full suspension for excellent outdoor performance

Wrangler Benefits

- 200-amp, Curtis microprocessor-based PMV controller system provides optimal levels of torque and power
- Two rugged independent in-line gear motors provide positive traction
- Variable speed up to 10 mph

Options

- Rear basket
- Safety flag
- Heavy-duty dust cover
- Weather breaker canopy
- Cup holder
- Rear basket

Silver Star Lifts & Ramps ★

Quad XL



This lift is compatible with the Hurricane® only.

Pride
Mobility Products Corp.

PRIDE
PERFORMANCE MOBILITY
— LIVE YOUR BEST —

Hurricane® & Wrangler®

STYLISH DESIGN &
PREMIUM PERFORMANCE.

PMV-PERSONAL MOBILITY VEHICLE

Pride's PMVs are a powerful blend of rugged, heavy-duty performance and refined luxury features. Whether you find yourself limited by a mobility impairment that threatens to interfere with your outdoor activities or you simply want a stylish, practical, and durable personal mobility vehicle that is easy to use, quiet, and comfortable, a Pride PMV is the perfect choice.



Pride®
Mobility Products Corp.

800-800-8586 (US) • 888-570-1113 (Canada) • www.pridemobility.com