



Q6 edge®

Q6
FEEL THE POWER®

Q6 EDGE® shown with Euro 2 Seating System and R-Net controller

.....

With features like standard 4-pole motors, Mid-Wheel 6® Drive Design and ATX Suspension, the Q6 Edge® is engineered to meet the performance needs of the most active user. The Q6 Edge accepts a complete range of seating and electronics options to deliver exceptional rehab capability at an excellent value.


LIFE BEYOND LIMITS™

www.quantumrehab.co.uk

Specifications

DRIVE WHEELS	35.56 cm (14") pneumatic or solid ³
LOW IMPACT OMNI-CASTERS	Front: 15.25 cm (6") solid Rear: 15.25 cm (6") solid
SUSPENSION	ATX Suspension (Active-Trac [®] with extra stability) with enhanced rear suspension
MAXIMUM SPEED ¹	Up to 4 mph (6.44 km/h)
BRAKES	Intelligent braking (electronic, regenerative disc brakes)
GROUND CLEARANCE	6.35 cm (2.5") at suspension arm 7.62 cm (3") at center of frame
TURNING RADIUS ²	50.8 cm (20")
OVERALL SIZE ²	Length: 90.17 cm (35.5") Width: 60.96 cm (24")
SEATING	Synergy [®] Seat Adjustable Width: 35.5 cm - 56 cm (14" -22") Adjustable Depth: 35.5 cm - 56 cm (14" -22") TRU-Comfort Plus Adjustable Width: 41 cm, 46 cm, 51 cm or 56 cm (16", 18" or 20") Adjustable Depth: 41 cm, 46 cm, 51 cm or 56 cm (16", 18" or 20") Euro 2 Adjustable Width: 41 cm, 46 cm or 51 cm (16", 18" or 20") Adjustable Depth: 41 cm, 46 cm or 51 cm (16", 18" or 20")
DRIVETRAIN	2-motor, 4-pole, Mid-Wheel 6 [®]
BATTERIES	Two 12 volt, deep cycle (NF-22) 55 AH
BATTERY CHARGER	8A, off-board
AVAILABLE ELECTRONICS	80 amp, PG R-Net LED controller
WEIGHT CAPACITY	21 Stone 6 lbs. (136 kg) maximum (user weight)
BASE WEIGHT	62 kg (138 lbs.)
BATTERY WEIGHT ⁴	17.24 kg (38 lbs.) (each)
WARRANTY	Lifetime limited warranty on frame; 1-year limited warranty on electronics; 18-month warranty on drive control system; 2-year limited warranty on drivetrain

Features

- Mid-Wheel 6[®] allows six wheels on the ground for maximum stability
- Compatible with TRU-Balance[®] Power Positioning
- Low impact OMNI-Casters (nylon, spherical-shaped casters) on front and rear prevent wheel hang-ups
- ATX Suspension (Active-Trac[®] with extra stability) incorporates front OMNI-Casters and semi-independent rear caster beam for enhanced performance over more varied terrain

Options

- TRU-Balance[®] Tilt and Recline
- 80 amp, R-Net LCD Colour controller
- Synergy[®] Seating System
- TRU-Comfort Plus Seating System
- TRU-Balance[®] Power Positioning System
- Cup holder
- Oxygen bottle holder
- Crutch holder
- Rear basket
- Leg rest options (Please see prescription form)
- Cantilever arm rests
- Synergy[®] centre mount adjustable foot platform
- Colour shrouds (for TRU-Comfort Plus only) ³
- Ottobock arm pads
- Gel arm pads




R-Net LED controller (standard)



R-Net Colour LCD controller (optional)

Shroud Colors

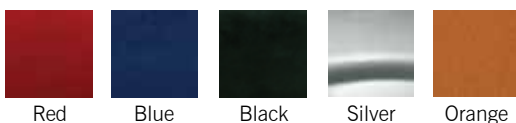

Red Blue Black Silver Orange

DEALER INFO HERE



32 Wedgwood Road, Bicester,
Oxfordshire OX26 4UL
(Tel) 01869 324 600
(Fax) 01869 323 070
www.quantumrehab.co.uk



FOOTNOTES:

- 1) Speed varies with user weight, terrain type, battery charge, battery condition and tyre pressure.
- 2) Due to manufacturing tolerances and continual product improvement, this specification can be subject to variance of + or - 3 %. Dimensions listed are for power base only. Overall measurements will vary based on seating and accessory selections.
- 3) Optional.
- 4) Battery weight may vary based on manufacturer.

The information contained herein is correct at the time of publication; due to the manufacturer's commitment to constant improvement and development, they reserve the right to alter specifications without prior notice. Theoretical range is calculated under test conditions in accordance with European standards. Actual range in normal use will depend on many factors, including the condition of the vehicle and its batteries, the weight of the driver, correct tyre pressure, ambient temperature, and the gradient and surface of the road or pavement.