

# Jazzy® 1420

The Jazzy® 1420 offers a balance of performance, styling and tight manoeuvrability. Designed especially for the bariatric client, the 1420 offers an impressive list of standard features including a weight capacity up to 227 kg (500 lbs.), 41 cm (16") knobby drive wheels, adjustable front anti-tips, Active-Trac® Suspension, and patented two-motor, mid-wheel drive to deliver superior outdoor performance without compromising indoor manoeuvrability.



Jazzy® 1420  
shown in Red

21 Healey Rd., Dandenong, Victoria 3175

ACN 088 609 661 ABN 66 088 609 661

(Tel) 03 9706 4611 • (Fax) 03 9706 4622 • [www.pridemobility.com](http://www.pridemobility.com)

**Pride**®  
Mobility Products  
Australia Pty. Ltd.

— LIVE YOUR BEST —

## SPECIFICATIONS

<b>Weight capacity</b>	227 kg (500 lbs.)
<b>Maximum speed</b>	Up to 7 km/h (4.5 mph)
<b>Ground clearance</b>	11 cm (4.25")
<b>Turning radius</b>	65 cm (25.75")
<b>Overall length</b>	117 cm (46")
<b>Overall width</b>	69 cm (27.25")
<b>Wheels</b>	Drive: 41 cm (16") pneumatic; mid-wheel
<b>Front anti-tips</b>	20 cm (8") solid; adjustable
<b>Rear castors</b>	23 cm (9") pneumatic; articulating
<b>Suspension</b>	Active-Trac®
<b>Drive train</b>	Two-motor, mid-wheel drive
<b>Braking system</b>	Intelligent braking (electronic, regenerative disc brakes)
<b>Standard electronics</b>	70 amp, PG Pilot Plus
<b>Specialty controls</b>	Yes
<b>Battery charger</b>	8 amp, Off-board (standard)
<b>Per-charge range</b>	Up to 32 km (20 miles)
<b>Battery requirements</b>	(2) 12 volt, deep cycle
<b>Battery size</b>	Group 24 (70 AH)
<b>Battery weight</b>	24 kg (53.5 lbs.) each
<b>Base weight</b>	64 kg (141.5 lbs.)
<b>Standard seat weight</b>	22.5 kg (49.6 lbs.) (euro-comfort)
<b>Std. seat dimensions</b>	46 cm x 46 cm (18" x 18")
<b>Warranty</b>	2-year limited

## STANDARD COLOUR CHOICES



Candy Apple  
Red




Viper Blue

## OPTIONAL COLOUR CHOICES



Purple

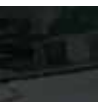

Black



Champagne

## FEATURES

- Independent rear and Active-Trac suspension for outstanding outdoor performance
- Adjustable front anti-tips
- Powerful, patented two motor mid-wheel drive design provides outstanding manoeuvrability
- Low centre of gravity design delivers a very smooth and stable ride
- 41 cm (16") knobby drive wheels and optimum weight distribution provide superior traction
- Speed range adjustable up to 9 km/h (6 mph) with high-speed motors
- Up to weight capacity of 227 kg (500 lbs.)
- Ergonomically enhanced seating

## OPTIONS

- Cup holder
- Rear basket
- Oxygen tank holder
- Cane/crutch holder
- Flat-free inserts
- Elevating leg rests
- Heel loops
- Adjustable-height armrests
- Multiple seating options
- Specialty controls
- Swing-away footrests

**Jazzy® 1420**  
shown in Red



21 Healey Rd., Dandenong, Victoria 3175

ACN 088 609 661 ABN 66 088 609 661

(Tel) 03 9706 4611 • (Fax) 03 9706 4622 • [www.pridemobility.com](http://www.pridemobility.com)

The information contained herein is correct at the time of publication; due to the manufacturer's commitment to constant improvement and development, they reserve the right to alter specifications without prior notice. Theoretical range is calculated under test conditions in accordance with European standards. Actual range in normal use will depend on many factors, including the condition of the vehicle and its batteries, the weight of the driver, correct tyre pressure, ambient temperature, and the gradient and surface of the road or pavement.

**PRIDE®**  
PERFORMANCE MOBILITY  
— LIVE YOUR BEST® —