

# Jazzy®

## 1450

The Jazzy® 1450 offers heavy-duty construction and performance components including upgraded motors, a reinforced frame, 90-amp PG VR2 controller and speeds up to 7.2 km/h (4.5 mph). This bariatric power chair features a 272 kg (600 lbs.) weight capacity, front-wheel drive for superb indoor maneuverability in the tightest of spaces and large front drive wheel for extra absorption over obstacles. The Jazzy 1450 combines superb climbing capabilities with front anti-tips for transitioning over grade changes and handling various terrains easily.



*Jazzy® 1450  
shown in Candy Apple Red  
with high-back seat*

**Pride**  
Mobility Products  
Australia Pty. Ltd.

(Tel) 03 8770 9600 • (Fax) 03 9703 2720 • [www.pridemobility.com.au](http://www.pridemobility.com.au)



## SPECIFICATIONS

<b>Weight capacity</b>	272 kg (600 lbs.)
<b>Maximum speed</b>	Up to 7.2 km/h (4.5 mph)
<b>Ground clearance</b>	3.81 cm (1.5")
<b>Turning radius</b>	67.3 cm (26.5")
<b>Overall length</b>	97.16 cm (38.25")
<b>Overall width</b>	74.3 cm (29.25")
<b>Drive wheels</b>	35.53 cm (14") solid
<b>Front wheels</b>	10.16 cm (4") solid anti-tips
<b>Rear wheels</b>	22.9 cm (9") castors
<b>Suspension</b>	Swing Arm
<b>Drivetrain</b>	Two-motor, front-wheel drive
<b>Braking system</b>	Intelligent braking (electronic, regenerative disc brakes)
<b>Standard electronics</b>	90 amp, PG VR2
<b>Battery charger</b>	Off-board, 8 amp
<b>Per-charge range</b>	Up to 19.62 km (12.19 miles)
<b>Battery requirements</b>	(2) 12 volt, Group 24
<b>Battery weight</b>	23.6 kg (52 lbs.) each
<b>Base weight</b>	88.5 kg (195 lbs.)
<b>Standard seat weight</b>	23.4 kg (51.5 lbs.) (high-back)
<b>Standard seat</b>	61 cm x 61 cm (24" x 24")
<b>Max. seat dimensions</b>	81 cm x 71 cm (32" x 28")

## FEATURES

- Black vinyl high-back seat with height-adjustable headrest
- Flat-free drive tires
- 152 cm (60") lap belt
- Heavy-duty, base mounted, flip-up height adjustable foot platform
- Dual-layer powder coated frame for increased durability

## OPTIONS & ACCESSORIES\*

- Pride® rear bag (includes crutch/ walking stick holder)
- Walker holder (must have Pride® rear bag to attach)
- Oxygen tank holder (must have Pride® rear bag to attach)
- Saddlebag
- Cup holder
- Mobile phone holder

\* Please note, accessory bracket may be needed for attaching rear bag.


## OPTIONAL SEATING

- Synergy® Seat
- Solid Seat Pan

## COLOUR CHOICES



Candy  
Apple Red

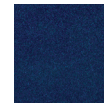

Viper  
Blue

**Pride**<sup>®</sup>  
Mobility Products  
Australia Pty. Ltd.

20-24 Apollo Drive Hallam, Victoria 3803  
ACN 088 609 661 ABN 66 088 609 661

(Tel) 03 8770 9600 • (Fax) 03 9703 2720 • [www.pridemobility.com.au](http://www.pridemobility.com.au)



The information contained herein is correct at the time of publication; due to the manufacturer's commitment to constant improvement and development, they reserve the right to alter specifications without prior notice. Theoretical range is calculated under test conditions in accordance with European standards. Actual range in normal use will depend on many factors, including the condition of the vehicle and its batteries, the weight of the driver, correct tyre pressure, ambient temperature, and the gradient and surface of the road or pavement.

