



Q6 edge **HD**

Get the Edge.

Q6 Edge® HD shown with
TRU-Balance® 3 HD Power
Positioning and Q-Logic 2
Drive Controls

With features like standard 4-pole motors, Mid-Wheel 6® Drive Design and ATX Suspension, the Q6 Edge® HD is engineered to meet the performance needs of the most active user. The Q6 Edge HD accepts a complete range of seating and electronics options to deliver exceptional rehab capability at an excellent value.



QUANTUM®

The Rehab Power Chair Solutions Company

www.quantumrehab.com

Specifications

DRIVE WHEELS	35.56 cm (14") solid (standard) or pneumatic (optional)
CASTOR WHEELS	Front: 15.24 cm (6") solid Rear: 15.24 cm (6") solid
SUSPENSION	ATX Suspension (Active-Trac® with extra stability) with enhanced rear suspension
MAXIMUM SPEED¹	Up to 7.6 km/h (4.7 mph)
BRAKES	Intelligent braking (electronic, regenerative disc brakes)
GROUND CLEARANCE²	7.62 cm (3") at base frame
TURNING RADIUS	50.2 cm (19.75")
SIZE²	Length: 90.4 cm (35.6") (without foot riggings) Width: 66.68 cm (26.25") (without foot riggings)
SEATING	Synergy® Seat Adjustable Width: 51 cm - 71 cm (20"-28") Adjustable Depth: 51 cm - 71 cm (20"-28") TRU-Balance® 2 HD Power Positioning Adjustable Width: 51 cm - 71 cm (20"-28") Adjustable Depth: 51 cm - 71 cm (20"-28") TRU-Balance® 3 HD Power Positioning Adjustable Width: 51 cm - 71 cm (20"-28") Adjustable Depth: 51 cm - 63.5 cm (20"-25") High-back Seat Width/Depth: 51 cm - 71 cm (20"-28") Medium-back Seat Width/Depth: 51 cm - 71 cm (20"-28")
DRIVETRAIN	Mid-Wheel 6®, 4-pole
BATTERIES	70 AH (Group 24)
RANGE PER CHARGE^{1,3}	Up to 25 km (15.5 miles)
BATTERY CHARGER	8A, off-board
AVAILABLE ELECTRONICS	100A Q-Logic Drive Control System
WEIGHT CAPACITY	204 kg (450 lbs.)
BASE WEIGHT	65 kg (143.5 lbs.) (without batteries)
BATTERY WEIGHT⁴	23.6 kg (52 lbs.) (each)
WARRANTY	Lifetime limited warranty on frame; 2-year warranty on electronics; 18-month warranty on drive motor

Features

- Mid-Wheel 6® allows six wheels on the ground for maximum stability
- Compatible with TRU-Balance® 3 HD Power Positioning
- Low impact OMNI-Castors (nylon, spherical-shaped castors) on front and rear prevent wheel hang-ups
- ATX Suspension (Active-Trac® with extra stability) incorporates front OMNI-Castors and semi-independent rear castor beam for enhanced performance over more varied terrain


Options

- TRU-Balance® 3 (450 lbs.) HD Tilt, HD Recline, and HD Tilt & Recline
- Headrest
- Occupied transit for TRU-Balance® 2 HD Tilt⁵
- Occupied transit for Synergy® Seating⁵

Shroud Colours



Candy Apple Red Viper Blue Onyx Black Sunburst Orange True Blue

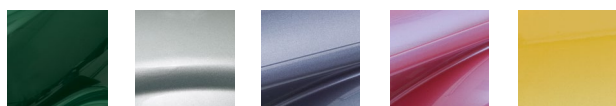

Emerald Green Silver Sapphire Blue Pearl Pink Orion Yellow

FOOTNOTES:

- 1) Range and speed vary with user weight, terrain type, battery charge, battery condition and tyre pressure.
- 2) Due to manufacturing tolerances and continual product improvement, this specification can be subject to variance of + or - 3%. Dimensions listed are for power base only. Overall measurements will vary based on seating and accessory selections.
- 3) Tested in accordance with ANSI/RESNA, WC Vol. 2, Section 4, standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Testing conducted at maximum weight capacities. This specification can be subject to a variance of +10%, -5%.
- 4) Battery weight may vary based on manufacturer.
- 5) Limited to 24" W x 24" D.

The information contained herein is correct at the time of publication; due to Pride's commitment to constant improvement and development, we reserve the right to alter specifications without prior notice.

Theoretical range is calculated under test conditions in accordance with European standards. Actual range in normal use will depend on many factors, including the condition of the vehicle and its batteries, the weight of the driver, correct tyre pressures, ambient temperature, and the gradient and surface of the road or pavement.

Pride®
Mobility Products
Australia Pty. Ltd.

— LIVE YOUR BEST™ —

AU_Q6 Edge HD_SS_7-09-14

20-24 Apollo Drive
Hallam, Victoria 3803
www.pridemobility.com.au

Tel: 03 8770 9600
Fax: 03 9703 2720


QUANTUM®

The Rehab Power Chair Solutions Company

