

Jazzy[®]

SELECT 6




The Jazzy Select[®] 6 incorporates sporty design and performance features for a truly unique ride experience. Pride's patented Active-Trac[®] Suspension combined with six wheels on the ground and an innovative articulating beam allow for superior performance indoors and out, with added stability over uneven terrain and ramp transitions.

ACTIVE-TRAC[®] SUSPENSION



**JAZZY ARMOR[®] AND
COLOR-THROUGH
SHROUD**

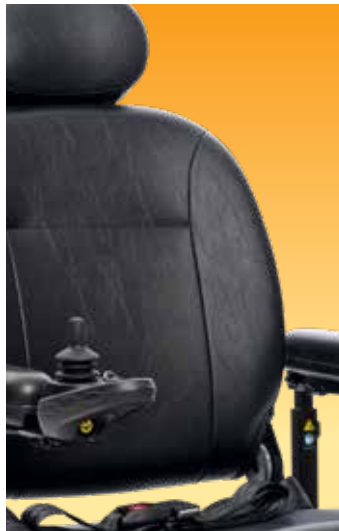

**COMFORT HIGH-BACK SEAT
WITH REPLACEABLE BACK
AND SEAT COVERS**

**NORTH AMERICA'S
FAVORITE POWER CHAIR**

Learn more online.



gojazzy.com

Specifications

Weight Capacity	300 lbs.
Maximum Speed	Up to 4 mph
Ground Clearance	1.75"
Turning Radius	21.875"
Overall Length	35.375"
Overall Width	23.5"
Drive Wheels	10"
Front Wheels	6" solid casters
Rear Wheels	6" solid casters
Suspension	Active-Trac® (ATX)
Drivetrain	Two-motor, in-line, mid-wheel drive
Braking System	Regenerative and electro-mechanical
Standard Electronics	40 amp, PG GC 3 controller
Battery Charger	Off-board, 3.5 amp
Per-Charge Range	Up to 15 miles
Battery Requirements	(2) 12 volt, U-1
Battery Weight	25 lbs. each
Base Weight	100 lbs.
Standard Seat Weight	43 lbs. (comfort, high-back)
Standard Seat	16" x 16-18", 18" x 18-20", 20" x 18-20"
Max. Seat Dimensions	20" x 20"
Warranty	5-year limited on frame; 13-month limited on batteries; 13-month limited on electronics; 13-month limited on drive motors

Colors



Jazzy Red



Jazzy Blue

Features and Benefits

- Jazzy® Armor® includes shroud and controller guards to protect against daily wear and tear
- Depth adjustable, limited recline comfort high back seat available
- Universal, fixed, in-line joystick mount
- Flip-up, angle-adjustable foot platform
- Black, flat-free, non-marking tires
- 60" lap belt

Seating Options*

- Depth adjustable, limited recline comfort high back seat with solid seat pan, available sizes: 16"x16"-18", 18"x18"-20", 20"x18"-20"
- Deluxe contour, limited recline high back available sizes: 18"x16"
- Deluxe contour, limited recline, high back with solid seat pan, available sizes: 16"x16"-20"x20"
- Power Elevating Seat (note seat does not swivel)
- Synergy® Seat

Options

- Bell shaped foot platform
- Elevating leg rests
- Heel loops
- Residual limb support
- Swing-away footrests
- Swing-away joystick
- VR2 6-key (must upgrade with PES option)

Accessories

- Cane/crutch holder
- Cell phone holder
- Cup holder
- Lap belt (50", 70", 80", 90")
- Oxygen tank holder
- Rear basket
- Saddle bag
- Walker holder
- Weather cover
- XLR USB charger

*Additional charges may apply.