

Litestream[®] XF LTD



SPECIFICATIONS

Product Weight:	27.5 lbs. ¹
Weight Capacity:	250 lbs.
Seat-to-Floor Heights:	12"-16.5" front 12"-16" rear
Seat Sizes:	Width: 15"-20" Depth: 14"-20" (Available in 1" increments)

COLORS

STANDARD



Black Abyss

OPTIONAL



Red Hot



Blue Sky

Beamsville, Ontario • 888-570-1113
Visit us on the web: www.quantumrehab.com

The Litestream[®] XF LTD ultralight manual wheelchair uses a unique, innovative, rigidizing double X folding frame that delivers the versatility and performance of a rigid frame while folding to a width of only 9" for easy transport. The Litestream XF LTD accepts a broad selection of accessories to meet a wide variety of needs and personal preferences while remaining an excellent overall value.

OPTIONS

REAR WHEELS:	20", 22", 24" & 26" Composite mag 20", 22", 24" & 26" Spinerger [®] wire spoke
TIRES:	20", 22", 24" & 26" Polyurethane tire 20", 22", 24" & 26" Everyday gray tire
HANDRIMS:	20", 22", 24" & 26" Aluminum 20", 22", 24" & 26" Vinyl-coated 20", 22", 24" & 26" Projection 20", 22" & 24" Chrome 22", 24" & 26" Natural-Fit (Standard and Supergrip)
AXLES:	Quick-release, quad-release
AXLE PLATES:	Adjustable axle plate (standard) Amputee mounting (optional)
CASTERS:	4", 5", 6" & 8" Black polyurethane
FORKS:	5" or 8" bolt-on or quick-release forks (angle-adjustable) Frog Legs Caster pin locks
WHEEL LOCKS:	Pull-to-lock Pull-to-lock with extensions Push-to-lock Push-to-lock with extensions
ARMRESTS:	Height-adjustable drop-in Dual-post, flip-back, height-adjustable arms Desk or full-length arm pads
BACK CANES:	8 degree with push handles Angle-adjustable, with push handles
BACKS:	Standard padded nylon sling with convenience pocket Tension-adjustable upholstery Synergy [®] Shape Back
FRONT RIGGINGS:	Aluminum, bi-directional swing-away leg rests, 60° & 70° Hemi Flip-up aluminum footplates (with heel loops) Composite comfort adjust footplates Angle-adjustable footplates Nylon heel loops High pivot elevating leg rests with standard or offset calf pad
ACCESSORIES:	Height-adjustable, rear anti-tip wheels Lap belt Tool kit Combat cane/crutch holder Unoccupied transport kit

¹) Standard configuration: 16" x 16", 19" backrest, sling upholstery, push-to-lock brakes, 22" composite wheels, flat-free tires & 6" casters (not including leg rests/armrests).

REAR WHEELS



20", 22", 24" & 26" composite mag wheels are exceptionally durable and low maintenance



20", 22", 24" & 26" high-performance Spinergy® wheels feature double-wall rims, with straight-pull spokes for reduced weight and optimal durability

CASTERS



4", 5", 6", & 8" composite casters feature precision bearings and long-life black polyurethane tires

TIRES



20", 22" & 24" Polyurethane Tire (Solid)



20", 22" & 24" Everyday Gray Tire (Pneumatic)

HANDRIMS



20", 22", 24" & 26" Aluminum Handrim



20", 22" & 24" Chrome Handrim (Only available on composite wheels)




20", 22", 24" & 26" Vinyl-coated Handrim



20", 22", 24" & 26" Projection Handrim (Available in vertical and 45°)

MISCELLANEOUS

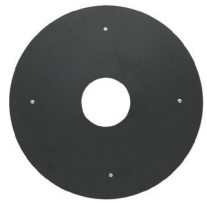

Spoke protectors

FORKS



8" tall fork



5" short fork

Available as high angle bolt-on and quick-release



Frog Legs suspension fork

WHEEL LOCKS



Pull-to-lock



Push-to-lock

ARMRESTS



Single-post height-adjustable armrests



Dual-post flip-back height adjustable arms

LEG RESTS



60° & 70° Hemi aluminum bi-directional swing-away leg rests



Flip-up aluminum footplates with heel loops



10-degree composite, comfort adjust footplates in 4.25", 5", and 6" widths



Angle-adjustable aluminum footplates in 5.5", 6.5" and 7.5" widths



Nylon heel loops

BACK OPTIONS



Tension-adjustable upholstery



Standard padded nylon sling with pocket



Synergy® Shape back

ACCESSORIES AND PERSONALIZATION



Combat, cane & crutch holder



Height-adjustable anti-tips



The Transit Securement Package allows convenient securement of the Litestream® XF LTD for safe transport. For additional information, please refer to the Transit Securement Package brochure.