

DEDICATED TO REVOLUTIONIZING THE POWER CHAIR EXPERIENCE



QUANTUM®

[quantumrehab.com](http://quantumrehab.com)



Front and rear LED fender lights for seeing and being seen (standard)



Drive wheel accents available in bold new colors (standard)



Conveniently located USB charger for recharging your gadgets on the go (standard)



# EDGE 3

The Edge 3 was designed to offer the most advanced power chair experience ever. Upgraded SRS (Smooth Ride Suspension) provides a smooth, comfortable ride and advanced stability. Optional iLevel® is available at 4.5 mph with up to 12" of power adjustable seat height. Standard features include front and rear LED fender lights, a USB charger, and drive wheel color accents. The Edge 3 is available in a total of 19 color options, including seven brand-new, bold colors.



**PROUDLY MADE IN**  
**AMERICA**   
OF FOREIGN AND DOMESTIC PARTS



Proudly made with care in America  
from foreign and domestic parts



# Q6 Series



Carbon Fiber

## Edge 3

- Featuring upgraded SRS (Smooth Ride Suspension)
- Available with iLevel® 12" power elevation at up to 4.5 mph (retrofitable on the Edge 3)
- Extra Stability Technology™ enhances safety while performing transfers, reach, and numerous other activities of daily living (available only with iLevel® option)
- Fender lights, USB charger, and drive wheel accents included
- 7 bold new colors increase the total available colors to 19
- Made in America of foreign and domestic parts





Yellow Blaze

## Q6 Edge® 2.0



- User weight capacity of 300 lbs.
- Mid-Wheel 6® allows six wheels on the ground for stability
- Standard 6 mph motors offer advanced performance, precise control and responsiveness via innovative single-stage gearbox
- iLevel® option available with up to 12" of power adjustable seat height at walking speed (3.5 mph)
- Available with Q-Logic 3 Advanced Drive Control System
- Up to 17 exciting color options
- Available with WC-19 crash-test approved occupied transit kit option on tilt & recline and iLevel®
- Available with Accu-Trac Advanced Tracking Technology
- Optional switch-operated front and rear LED fender lights (standard on iLevel®)
- Optional 8" casters



Lime Light

## Q6 Edge® Z



- User weight capacity of 300 lbs.
- Mid-Wheel 6® allows six wheels on the ground for stability
- 8 mph motor package option available
- iLevel® option available with up to 12" of power adjustable seat height at walking speed (3.5 mph)
- 13 exciting color options
- Switch-operated front and rear LED fender lights (standard)
- Available with Accu-Trac Advanced Tracking Technology



Electric Blue

## Q4

- User weight capacity of 300 lbs.
- Mid-Wheel 6® allows six wheels on the ground for stability
- Standard 6 mph motors offer advanced performance, precise control and responsiveness
- Bold aesthetics
- Q-Logic 2 NE+ electronics standard
- High-back seat
- ATX Suspension (Active-Trac® with extra stability) incorporates front OMNI-Casters and semi-independent rear caster beam for enhanced performance over more varied terrain
- TRU-Balance® Flex seating



## J6

- Mid-Wheel 6® allows six wheels on the ground for maximum stability
- Compatible with TRU-Balance® 2 Tilt
- Side-mounted, easily accessible freewheel levers
- ATX Suspension (Active-Trac® with extra stability) incorporates enhanced performance over more varied terrain
- Easy front access to batteries



# Front-wheel Drive

## 4Front®

The 4Front® is a quiet, more responsive front-wheel drive power chair that features automotive grade suspension with unprecedented comfort and ride quality. The 4Front's CASE (Caster Angle Sensor Encoder) technology and SRS (Smooth Ride Suspension) provide stability for smooth obstacle transitions and unmatched driving performance. The 4Front features standard speeds of up to 6 mph in the lowered position, and speeds up to 3.2 mph while elevated at 10" with our optional safe seat elevation. Standard LED fender lights allow you to see and be seen.

### FEATURES

- CASE (Caster Angle Sensor Encoder) technology
- Unmatched driving performance and ride quality
- Advanced stability for smooth obstacle transitions
- Superior, automotive-grade SRS (Smooth Ride Suspension)
- Available with optional 10" of power adjustable safe seat elevation at speeds of up to 3.2 mph
- Q-Logic 3 EX Advanced Drive Controls
- LED front and rear fender lights
- Quiet ride technology
- Bold aesthetics with 17 color options
- Available with WC-19 crash-test approved, occupied transit option with all configurations







# iLevel® Power Seat Elevation



12" LIFT  
AT UP TO  
4.5 MPH

Maximum lift of 10" on Q6 Edge® HD

12" of power adjustable seat height while  
at walking speed (10" on Q6 Edge® HD)

iLevel's patented Extra Stability Technology™  
enhances safety for transfers, reach, and  
numerous other activities of daily living

Up to 19 exciting color options

Retrofittable on the Edge 3 and Q6 Edge® 2.0

300 lbs. weight capacity on Edge 3,  
Q6 Edge® 2.0, and Q6 Edge® Z

400 lbs. weight capacity  
on Q6 Edge® HD

Extra Stability Technology™ for  
a smooth and stable ride

Front and rear LED fender lights  
come standard with iLevel®



iLevel® Power Chairs feature a revolutionary new power chair design for a better experience. iLevel improves daily functionality and quality of life by allowing operation of the power chair with the seat elevated while at walking speed (up to 4.5 mph).



# QUANTUM<sup>®</sup> KOZMO<sup>™</sup>

and Kozmo<sup>™</sup> Portable

## FEATURES

- Max. speed up to 3.6 mph
- 200 lbs. weight capacity
- One-hand feather-touch disassembly available on Kozmo<sup>™</sup> Portable
- Choice of 6 colors
- 18 amp batteries
- Compact turning radius of 25.45"
- Pediatric, height-adjustable footplate
- 38" lap belt



The Kozmo<sup>™</sup> Portable can be disassembled quickly for transport or storage.



## Color Options



Ruby Red



Sapphire Blue



Rose Quartz



Citrine Yellow



Amber Orange



Pearl White

## Pediatric

Quantum<sup>®</sup> is proud to offer power chairs with features such as power seat elevation, SRS (Smooth Ride Suspension), LED lighting, Bluetooth<sup>®</sup>, and color options that not only appeal to children but give parents peace of mind. With seating sizes as small as 10" x 10", Quantum pediatric power chairs are designed to grow with your child, offering comfort and reliability for years to come.



## FEATURES\*

- Available with power seat elevation technology
- Superb driving performance and ride quality
- Maximum speeds up to 6 mph
- LED front and rear fender lights to see and be seen
- More than a dozen kid-friendly color options
- WC-19 crash-test approved, occupied transit option (Q6 Edge<sup>®</sup> 2.0 and 4Front<sup>®</sup> only)
- Q-Logic 3 environmental controls featuring Bluetooth<sup>®</sup>
- Made in America of foreign and domestic parts

## Edge 3

Seat sizes as small as 12" x 12"

Available with 12" iLevel<sup>®</sup> (4.5 mph)

19 color options



## Q6 Edge<sup>®</sup> 2.0

Seat sizes as small as 12" x 12"

Available with 12" iLevel<sup>®</sup> (3.5 mph)

17 color options



## 4Front<sup>®</sup>

Seat sizes as small as 12" x 12"

Optional 10" elevation (3.2 mph)

17 color options



\* Some features are not available on all models. Ask your Authorized Quantum Dealer or visit [quantumrehab.com](http://quantumrehab.com) for details.



## Rear-wheel Drive



### Quantum® Rival®

- 12 exciting color options
- Compact size for outstanding maneuvering in tight spaces
- Maximum speed up to 8 mph with high-speed motor package
- High-performance, high efficiency motor package
- Accepts a wide range of seating and positioning systems
- Two suspension packages
- Compatible with TRU-Balance® 3 Power Positioning Systems
- Sliding battery tray for easy access

## Synergy® Seating

The Synergy® Seating System offers a clean, streamlined appearance and unparalleled ease of adjustment. The seating system's unique design allows for height, width, depth and back angle adjustments to be made quickly and easily in the field. The Synergy Seat interfaces with most types of seat and back cushions, and is designed to be the most user-friendly system available today.

### FEATURES

- Easily adjustable and versatile design that can grow with the client
- Easy customization
- Several width/depth adjustment packages ranging from 10" to 32"
- Accepts a wide variety of after-market seats and backs
- Heavy Duty option provides a weight capacity of 600 lbs. and seating up to 32" x 28" (model specific)







Emerald Green

## Quantum® 1450

- 600 lbs. weight capacity
- Maximum speed up to 4.5 mph
- High-performance, high-efficiency motor package
- Compatible with TRU-Balance® Power Positioning and a wide range of other seating and positioning systems



Midnight Blue

## Q6 Edge® HD

- Compatible with TRU-Balance® 3 HD Power Positioning
- ATX Suspension (Active-Trac® with extra stability) incorporates front OMNI-Casters and semi-independent rear caster beam for enhanced performance over more varied terrain
- 450 lbs. weight capacity (400 lbs. with iLevel® option)
- iLevel® option available with up to 10" of power adjustable seat height at walking speed (3.5 mph)
- Mid-Wheel 6" for maximum stability
- Optional switch-operated front and rear LED fender lights (standard on iLevel®)



## TRU-Balance® Flex Seating

TRU-Balance® Flex Seating is the no-charge static seat option available on the Edge 3, Q6 Edge® 2.0, and the Rival® power chairs.

### FEATURES

- Easily adjustable, this versatile design can grow with the user
- Same ease of width and depth adjustment as TRU-Balance® 3
- Height adjustable single post arm with clothing guard standard
- Accepts all rail-mounted TRU-Balance® 3 accessories
- Seat sizes 15"-22" wide and 15"-20" deep
- Back canes only

# TRU-Balance® 3



## Power Recline

- 85° - 175° of recline
- 300 lbs. weight capacity
- Does not increase seat-to-floor height
- Standard shear reduction and cantilever armrests
- Available in up to 19 back shroud colors
- Available with TRU-Balance® 3 Power Tilt and TRU-Balance® 3 Power Adjustable Seat Lift



## Power Tilt

- 50° of tilt
- Lower standard seat-to-floor height
- Unmatched adjustment and durability
- Maintains center of gravity shift
- Available in up to 19 back shroud colors



## Power Adjustable Seat Lift with iLevel®

- iLevel® technology enables the unit to drive up to 4.5 mph while fully elevated
- Offers 12" of power adjustable seat lift (10" maximum lift on Q6 Edge® HD)
- Available in up to 19 back shroud colors
- Available on the Edge 3, Q6 Edge® 2.0, Q6 Edge® Z, and Q6 Edge® HD



## Power Articulating Foot Platform

- Broad articulation range for enhanced adjustment
- 450 lbs. weight capacity
- Unmatched adjustment, durability and strength
- Multiple calf pads, calf panel and footplate options



## Heavy Duty Power Tilt & Recline

- 450 lbs. weight capacity
- 50° of tilt
- 85° - 175° of recline
- Available in up to 13 back shroud colors



## Complimentary Cup Holder

**Please Note:** Some options may not be available on all power bases. Ask Quantum Customer Service for details.





## TRU-Balance® 3 Power Positioning Systems

- Make fine adjustments while client is in the chair
- Order your seat width independent of your back width at no additional charge
- Internally or externally rotate cantilever armrests for optimal positioning
- Easily retrofit power tilt and power recline
- iLevel® is retrofittable on the Edge 3 and Q6 Edge® 2.0
- Easy back height and depth adjustment
- Easy center of gravity adjustment
- Complimentary cup holder
- Premium accessories available
- Optional attendant control
- Tool kit included for quick adjustments



# Cushions and Backs



## Cushion Color Guide

To make your cushion selection process a little easier, Stealth Products® Cushions are all marked with a color-coded PVC patch. The color of the patch identifies the function and HCPC code of the cushion.



Stealth Products®  
Simplicity G

- Contoured, high-density molded foam base
- Coccyx and sacral relief
- Moisture-wicking and easy-to-clean Stealth-Tek incontinence cover
- Sizes available: 16"W x 16"D - 24"W x 18"D



Stealth Products®  
Spectrum® Foam SPP

- Lateral and pre-ischial shelf contoured, high-density molded foam base
- Multi-layer, viscoelastic foam insert under at risk bony areas
- CoolCore® technology cover standard
- Urethane liner for incontinence
- Optional four-way stretch poly liner
- Sizes available: 16"W x 16"D - 32"W x 28"D
- Optional Stealth-Tek or mesh cover



Stealth Products®  
Spectrum® Gel SPP

- Lateral and pre-ischial shelf contoured, high-density molded foam base
- Twin-cell, lightweight fluid insert under at risk bony areas
- CoolCore® technology cover standard
- Urethane liner for incontinence
- Optional four-way stretch poly liner
- Sizes available: 14"W x 16"D - 24"W x 24"D
- Optional Stealth-Tek or mesh cover



Stealth Products®  
Solution® SPP

- Lateral and pre-ischial shelf contoured, high-density molded foam base
- Twin-cell, aqueous fluid insert under at risk bony areas
- CoolCore® technology cover standard
- Urethane liner for incontinence
- Optional four-way stretch poly liner
- Sizes available: 14"W x 16"D - 24"W x 24"D
- Optional Stealth-Tek or mesh cover



Active Contour  
ADI Q-Back



Extra Deep Contour  
ADI Q-Back



Deep Contour  
ADI Q-Back

## ADI Q-Back

The ADI Q-Back from Stealth Products®, Inc. is available on TRU-Balance® 3 Seating Systems and is offered in three configurations:

- Active Contour ADI Q-Back
- Deep Contour ADI Q-Back
- Extra Deep Contour ADI Q-Back

### Available sizes for tilt and static seating:

Heights: 19"-22"  
Widths: 16", 18", 20"


### Available sizes for power recline seating:

Heights: 22"-25"  
Widths: 16", 18", 20"

*Note: Additional back sizes are available with Special Order.*



# Cushions and Backs

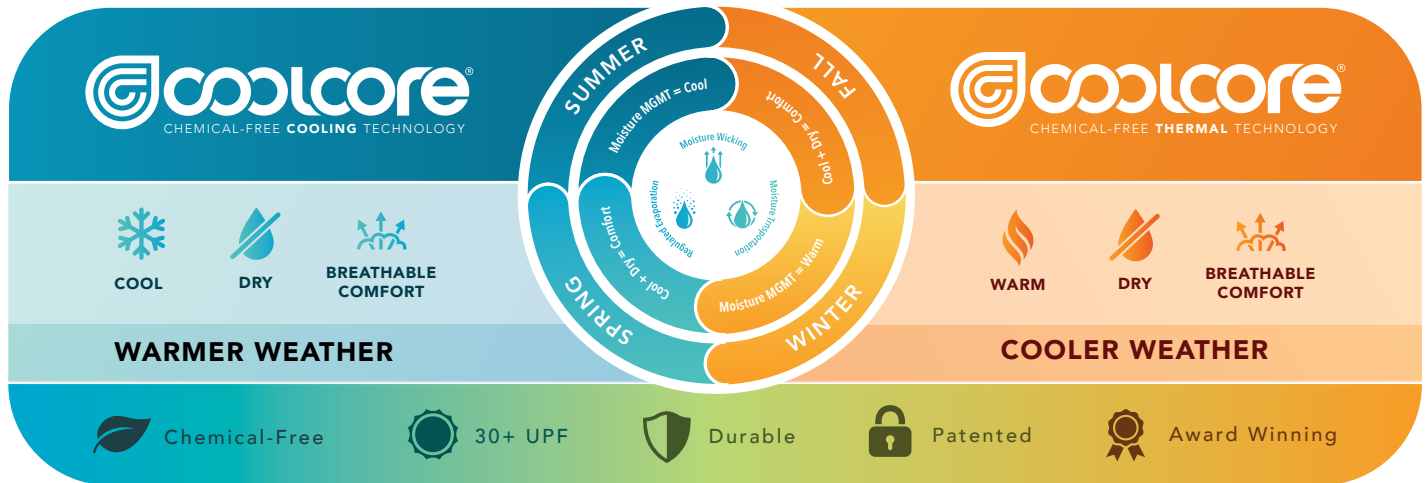

**Note:** CoolCore® is not available on the ADI Q-Back.



## SPP

• Lateral shelf contoured, high-density specially formulated molded foam base  
• Gel insert under at coccyx  
• Standard Visco cover standard for added comfort  
• Sizes available: 16"W x 16"D - 24"W x 24"D  
• Custom sizes available  
• Optional mesh cover



## Stealth Products® TRU-Comfort 2 SPP

- Deep contoured, high-density molded foam base
- CoolCore® technology cover standard
- Easy-pull rubber loop zipper handles
- Convenient side pouches on both sides of the cushion
- Sizes available: 16"W x 14"D - 22"W x 24"D
- Custom sizes available
- Optional four-way stretch poly liner




## Stealth Products® Essence SPP

- Available in Basic-Contour, Low-Contour and Medium-Contour configurations
- Lateral and pre-ischial shelf contoured, high-density specially formulated molded foam base
- Standard Visco topper for added comfort
- Coccyx and sacral relief
- CoolCore® technology cover standard
- 4-Way Stretch Poly liner optional
- Sizes available - 16"W x 16"D - 20"W x 20"D
- Custom sizes available
- Optional Mesh outer cover



## Stealth Products® Essence Plus SPP

- Standard adjustable high profile wedge kit
- Lateral and pre-ischial shelf contoured, high-density specially formulated molded foam base
- Standard Visco topper for added comfort
- Coccyx and sacral relief
- CoolCore® technology cover standard
- 4-Way Stretch Poly liner optional
- Sizes available - 16"W x 16"D - 20"W x 20"D
- Custom sizes available
- Adjustable contour wedge package provides better positioning customization
- Topper options include, gel or HR foam
- Optional Mesh outer cover

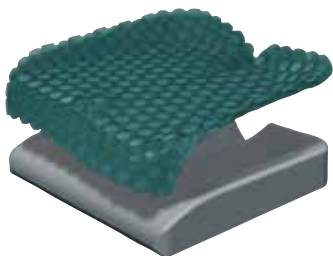

## Synergy® Solution® 1

- Deep contoured, high-density molded foam base
- Coccyx and sacral relief
- Soft gel top layer for skin protection
- Four-way stretch, fluid-resistant upper cover



## TRU-Comfort 2 Back Cushion

- Adjustable positioning and pressure management
- Compatible with both TRU-Balance® 3 Power Positioning Systems and the Static Seat Frame System
- Four-way stretch, recovery, fluid-resistant, breathable cover

# Electronics

Featuring  
 **Bluetooth®**  
(standard)



Q-Logic 3 EX  
Joystick Controller



Q-Logic 3 EX  
Enhanced Display

## Q-Logic 3 Drive Control System

The Q-Logic 3 Advanced Drive Control System from Quantum Rehab® is like no other power wheelchair electronics system. Designed from clinician and consumer feedback, Q-Logic 3 is user-friendly and highly customizable. This adaptable, expandable control system can increase and maximize independence.



Home Screen



Drive Screen

### FEATURES

- Wireless Bluetooth® programming through Windows devices
- Bluetooth access to computer, tablet, and phone functions
- iAccess programmable push button and toggle options for customized seating access
- Simple firmware updates
- Capable of utilizing multiple input devices
- Preset drive settings
- Preset color background options
- Clinic Mode:
  - Quick setup of multiple input devices during client evaluations
  - Automatically assigned drive modes to input devices

## Q-Logic 3 EX Enhanced Display

The optional Q-Logic 3 EX Enhanced Display serves as a separate screen used with specialty drive controls to display system data. With the Enhanced Display's built-in environmental controls, Q-Logic 3 can operate a wide range of appliances that use infrared (IR) remote controls, such as televisions, and serves as a Bluetooth® computer mouse emulator.



## Optional Modules



Q-Logic 3 SAJ



Q-Logic 3 iAccess



Q-Logic 3 Specialty Controls  
Interface Module





Jesse Cuellar, Artist

## Non-Expandable



Q-Logic 2  
NE and NE+



PG VR2 (4-key)



PG VR2 (6-key)



PG VR2 (9-key)

## Expandable (Q-Logic 3)



Q-Logic 3 EX  
Joystick



Q-Logic 3 EX  
Enhanced  
Display



Q-Logic 3 EX  
SAJ



Q-Logic 3  
Specialty  
Controls  
Interface  
Module

## Expandable (Q-Logic 2)



Q-Logic 2 EX  
Joystick



Q-Logic 2 EX  
Enhanced  
Display



Q-Logic 2  
EX Stand Alone  
Joystick /  
Attendant  
Joystick



Q-Logic 2  
Specialty  
Controls  
Interface  
Module

## Multi-Switch Systems



Stealth® i-Drive®  
Tri-Array Series,  
3-Switch



Stealth® Ultra  
3-Switch Head  
Array

## Proportional Drive Controls



Mo-Vis™  
All-Round  
Joystick



Mo-Vis™  
All-Round Light  
Joystick



Mo-Vis™  
Multi Joystick



Mo-Vis™  
Micro Joystick



Stealth®  
Precision Mini  
Proportional  
Joystick



Stealth®  
Mushroom  
Joystick

## Single Switch



Stealth®  
Egg Switch



Mo-Vis™  
Twister Basic  
Button



Mo-Vis™  
Twister Satellite  
Button

## Sip-N-Puff



































Stealth®  
Sip-N-Puff

## Misc. Options



Joystick  
Handles

# Color Options

		Edge 3	Edge 2.0	4Front®	Edge Z	Edge HD	Rival®	Q4	J6	Q1450
Sapphire Blue										•
Deep Purple		•	•	•	•	•	•			
Sugar Plum		•								
Midnight Blue		•	•	•	•	•	•			
Blue Jean		•								
Viper Blue									•	•
Electric Blue		•	•	•	•	•	•	•		
True Blue										•
Emerald Green										•
Military Green		•	•	•	•	•				
Night Train		•	•	•						
Onyx Black										•
Back in Black		•	•	•	•	•	•	•		
Grey Street		•	•	•	•	•	•			
Silver										•
Silver Bullet			•	•	•	•	•			
White Satin		•	•	•	•	•	•			
Green Machine		•								
Lime Light			•	•	•	•	•			
Yellow Blaze			•	•	•	•	•			
Lemon Crush		•								
Orion Yellow										•
Sunburst Orange										•
Orange Crush			•	•	•	•	•			
Evolution Orange		•								
Candy Apple Red									•	•
Red Corvette		•	•	•	•	•	•	•		
Cherry Bomb		•								
Raspberry Beret		•								
Pearl Pink										•
Pretty in Pink			•	•	•	•	•			
Carbon Fiber		•								
Digital Camo		•								



# Service & Support



stealthproducts.com

Stealth Products® launched in 1999, in response to a need for versatile, high-quality head and neck positioning systems, as well as switch site location hardware. Their products include: head and body support/positioning, joystick mounts, trays, seating, and pediatric mobility products.

## Education (888-570-1113)

Contact your Pride representative or Pride Technical Service for additional information on training available.

## Sales Support (888-570-1113)

Assistance with design, ordering, pricing and quality assurance.

## Custom Group (888-570-1113)

Providing innovative custom solutions and specialized product modifications – your “off the order form” solution team. Contact your Quantum representative or Pride Technical Service for assistance and additional information

## Technical Support (888-570-1113)

Assistance with troubleshooting and parts order placement.

## Pride Demo Program


Now available online. Contact customer service or your Pride representative for additional information

## Quick Ship Programs

Contact your Pride representative or your customer service representative for additional information.

## Set to Spec Service

Set to Spec Service from Quantum Rehab® is designed to assist you in delivering Q6 Power Chairs to your clients. We will pre-set, to your specifications, seat-to-floor heights, seating center of gravity, armrest heights, front riggings, thoracic and/or thigh laterals.

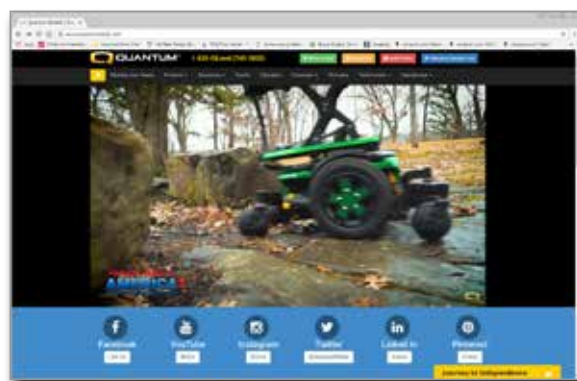

## quantumrehab.com

The ultimate resource for everything you need to know about Quantum® products.



## Clinical Resources

A special section of [quantumrehab.com](http://quantumrehab.com) for clinicians, featuring order forms, clinical references, and reimbursement information.



## Follow us...

Stay up to date on the latest Quantum® events and product releases through social media:



# Litestream® Manual Wheelchairs



The Litestream® XF ultralight manual wheelchair uses a unique, innovative, rigidizing double X folding frame that delivers the versatility and performance of a rigid frame while folding to a width of only 9" for easy transport. The Litestream XF accepts a broad selection of accessories to meet a wide variety of needs and personal preferences while remaining an excellent overall value. The Litestream XF is also available in an HD configuration that weighs just 35.5 lbs. and features a cane-to-cane rigidizing bar for added stability.



**Litestream® XF**  
with HD configuration



**Litestream® XF**  
with Ultra-hemi configuration



**Litestream® XF**  
with standard configuration



The Litestream® XF LTD ultralight manual wheelchair uses a unique, innovative, rigidizing double X folding frame that delivers the versatility and performance of a rigid frame while folding to a width of only 9" for easy transport. The Litestream XF LTD accepts a broad selection of accessories to meet a wide variety of needs and personal preferences while remaining an excellent overall value.



**Litestream® XF LTD**




## SPECIFICATIONS

	<i>Product Weight:</i>	<i>Weight Capacity:</i>	<i>Seat-to-Floor Heights:</i>	<i>Seat Sizes:</i>
<b>Litestream® XF</b>	26.5 lbs. <sup>1</sup>	250 lbs.	15"-19.5" (front) 14.5"-19" (rear)	Width: 14"-21" Depth: 14"-20"
<b>Litestream® XF HD</b>	35.5 lbs. <sup>3</sup>	350 lbs.	15"-19.5" (front) 14.5"-19" (rear)	Width: 18"-22" Depth: 18"-20"
<b>Litestream® XF Ultra Hemi</b>	26.5 lbs. <sup>1</sup>	250 lbs.	12.5"-16.5" (front) 12.5"-19" (rear)	Width: 14"-21" Depth: 14"-20"
<b>Litestream® XF LTD</b>	27.5 lbs. <sup>2</sup>	250 lbs.	12.5"-18.5" (front) 12"-17" (rear)	Width: 14"-21" Depth: 14"-20"


## COLORS

### LITESTREAM® XF, HD AND ULTRA HEMI

#### METALLIC

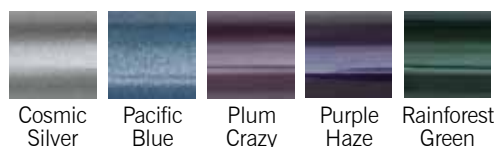

#### WATER TRANSFER PRINTING

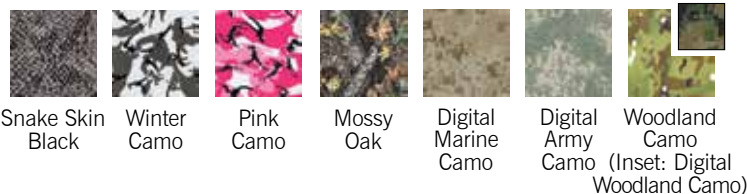

#### NON-METALLIC




### LITESTREAM® XF LTD

#### STANDARD



#### OPTIONAL

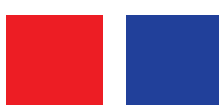

#### FOOTNOTES:

- <sup>1</sup> Standard configuration: 16" x 16", 19" back rest, sling upholstery, push-to-lock brakes, 24" spoke wheels, 1-3/8" pneumatic tires, 6" casters (not including leg rests).
- <sup>2</sup> Standard configuration: 16" x 16", 19" backrest, sling upholstery, push-to-lock brakes, 22" composite wheels, flat-free tires & 6" casters (not including leg rests/armrests).
- <sup>3</sup> Standard configuration: 18" x 18", 19" steel back canes, sling upholstery, push-to-lock brakes, 24" spoke wheels, solid tires, 6" casters (not including leg rests),  
HD cross frame and back cane rigidizing bar.



ALSO AVAILABLE:  
KOZMO™ PORTABLE



AVAILABLE  
WITH 12" iLEVEL®



AVAILABLE  
WITH 12" iLEVEL®



	KOZMO™	J6	Q4	Q6 EDGE® 2.0	EDGE 3
DRIVE CONFIGURATION	Rear-wheel Drive	Mid-wheel Drive	Mid-wheel Drive	Mid-wheel Drive	Mid-wheel Drive
DRIVE WHEELS	9" Solid	10" Solid Optional: 10" Pneumatic	14" Solid Optional: 14" Pneumatic	14" Solid Optional: 14" Pneumatic; 14" Knobby	14" Solid Optional: 14" Pneumatic
CASTER WHEELS	Front: 6" Solid	Front: 5" Solid Rear: 6" Solid	Front: 6" Solid Rear: 6" Solid	Front: 6" Solid (Optional: 8") Rear: 6" Solid (Optional: 8")	Front: 6" Solid Rear: 6" Solid
ANTI-TIP WHEELS	3" Solid	N/A	N/A	N/A	N/A
MAXIMUM SPEED <sup>4</sup>	Up to 3.6 mph	Up to 4 mph	Up to 6 mph	Up to 6 mph	Up to 6 mph (4)
GROUND CLEARANCE	1.6" at anti-tip bracket	2.5" (frame)	2.5" (at caster pivot)	2.7" (battery tray)	2.7" (battery tray)
TURNING RADIUS <sup>5</sup>	25.45" with foot platform	22"	20"	6" casters: 20.5" 8" casters: 23.25"	20.5"
OVERALL LENGTH <sup>11</sup>	33.09" with foot platform	34.74"	35.5"	6" casters: 32.23" 8" casters: 40.48"	32.23" without foot platform
BASE WIDTH	22.6"	23.19"	24"	24.4"	24.4"
<b>Seating Sizes:</b>					
SYNERGY®	N/A	W: 10-20" D: 10-20"	N/A	W: 10-22" D: 10-22"	N/A
TRU-BALANCE® 3 STATIC	N/A	W: 10-20" D: 10-20"	N/A	W: 10-22" D: 10-22"	W: 10-22" D: 10-22"
TRU-BALANCE® 3 POWER TILT	N/A	W: 14-20" D: 14-20"	N/A	W: 12-22" D: 12-22"	W: 12-22" D: 12-22"
TRU-BALANCE® 3 LIFT & TILT	N/A	N/A	N/A	W: 12-22" D: 12-22"	W: 12-22" D: 12-22"
TRU-BALANCE® 3 HD STATIC	N/A	N/A	N/A	N/A	N/A
TRU-BALANCE® 3 HD POWER TILT	N/A	N/A	N/A	N/A	N/A
TRU-BALANCE® 3 HD LIFT & TILT	N/A	N/A	N/A	N/A	N/A
TRU-BALANCE® 2 HD STATIC	N/A	N/A	N/A	N/A	N/A
TRU-BALANCE® 2 POWER TILT	N/A	N/A	W: 16-20" D: 16-20"	N/A	N/A
TRU-BALANCE® 2 HD POWER TILT	N/A	N/A	N/A	N/A	N/A
TRU-BALANCE® FLEX SEAT	N/A	N/A	W: 15-22" D: 15-20"	W: 15-22" D: 15-20"	W: 15-22" D: 15-20"
KOZMO™ PEDIATRIC STATIC SEAT	W: 12-16" D: 12-16"	N/A	N/A	N/A	N/A
<b>Seat-to-Floor Heights <sup>16,17</sup>:</b>					
SYNERGY® / STATIC <sup>8</sup>	N/A	17.75"	N/A	16.5"-18.5"	16.5" - 18.5"
POWER TILT <sup>8</sup>	N/A	17.75"	17.3"-19.3"	17.5"-19.5"	17.5" - 19.5"
LIFT & TILT <sup>8</sup>	N/A	N/A	N/A	17.5"-19.5"	17.5" - 19.5"
TRU-BALANCE® FLEX SEAT <sup>8</sup>	N/A	N/A	16.5"-18.25"	16.5"-18.25"	16.5"-18.25"
KOZMO™ PEDIATRIC STATIC SEAT <sup>8</sup>	15.75"-16.75"	N/A	N/A	N/A	N/A
MANUAL RECLINE	No	Yes	Yes	Yes	Yes
BATTERY SIZE <sup>6</sup>	(2) 18A, 12 V, deep cycle	U-1	NF-22	NF-22	NF-22 (2)
BATTERY WEIGHT <sup>14</sup>	31 lbs.	23.4 lbs.	39 lbs.	38 lbs.	38 lbs. each
AVAILABLE ELECTRONICS <sup>12</sup>	40A, PG GC3 (non-programmable)	70A Q-Logic 2 NE 70A Q-Logic 2 NE+ 70A Q-Logic 2 EX	120A Q-Logic 2 EX 75A Q-Logic 2 NE+	120A Q-Logic 3 EX 75A Q-Logic 2 NE 75A Q-Logic 2 NE+	120A Q-Logic 3 EX 75A Q-Logic 2 NE 75A Q-Logic 2 NE+
BATTERY CHARGER	3.5A Off-board	5A Off-board	8A Off-board	8A Off-board	8A Off-board
MOTOR PACKAGES	N/A	2-Pole 4 mph	4-Pole Mid-Wheel 6®	2-Pole Mid-Wheel 6®	2-Pole Mid-Wheel 6®
WEIGHT CAPACITY <sup>7</sup>	200 lbs.	300 lbs. <sup>9</sup>	300 lbs.	300 lbs.	300 lbs.
BASE WEIGHT	67 lbs.	100.2 lbs.	137 lbs.	138.6 lbs. (6" casters) 148 lbs. (8" casters)	138.6 lbs.

1) Standard  
2) Optional  
3) Requires HD Suspension  
4) Speed varies with user weight, terrain type, battery charge, battery condition and tire pressure.  
5) Without front riggings  
6) Battery capacity may change when adding power positioning systems.

7) Weight capacity is dictated by the selection of seating.  
8) All measurements +/- 1/8" with 18" x 18" seat, fixed towers, 300 lbs. weight capacity and solid tires (add 3/8" with pneumatic tires). Depths may increase seat-to-floor height. Contact Customer Service for details. Listed measurements are in the lowest position; does not include elevating system unless specified.  
9) 250 lbs. weight capacity with TRU-Balance® Systems  
10) A minimum seating size of 20" x 20" is required with iLevel® on the Q6 Edge® HD; maximum seating size is 24" x 25"

11) Without foot riggings  
12) All Quantum units can be ordered with specialty controls  
13) iLevel® not available with 8 mph motor package  
14) Battery weight may vary +/- 3 lbs. based on manufacturer.  
15) Weight capacity is reduced to 400 lbs. with iLevel® power adjustable seat height

16) Seat-to-floor height  
Due to the 4Front® design, seat-to-floor heights are approximately 16.5" seat-to-floor or iLevel®. 19.5" static seat or power



AVAILABLE WITH  
10" HD iLEVEL® AT 4.5 MPH



AVAILABLE WITH 12" iLEVEL®



AVAILABLE WITH 10" ELEVATION



AVAILABLE WITH 10" HD iLEVEL®



	Q6 EDGE® Z	4FRONT®	QUANTUM® RIVAL®	Q6 EDGE® HD	QUANTUM® 1450
Drive	Mid-wheel Drive	Front-wheel Drive	Rear-wheel Drive	Mid-wheel Drive	Front-wheel Drive
Pneumatic	14" Run Flat <b>Optional:</b> 14" Pneumatic; 14" Knobby	14" Flat Free <b>Optional:</b> 14" Knobby	14" Solid <b>Optional:</b> 14" Pneumatic; 14" Knobby	14" Run Flat <b>Optional:</b> 14" Pneumatic; 14" Knobby	14" Solid
	<b>Front:</b> 6" Solid <b>Rear:</b> 6" Solid	<b>Rear:</b> 9" Flat Free	<b>Front:</b> 9" Pneumatic	<b>Front:</b> 6" Solid <b>Rear:</b> 6" Solid	9" Solid
	N/A	<b>Front:</b> 4"	N/A	N/A	4" Solid
4.5 mph w/iLevel)	Up to 6 mph (Optional: 8 mph) <sup>13</sup>	Up to 6 mph	Up to 6 mph (Optional: 8 mph)	Up to 4.7 mph	Up to 4.5 mph
Ray)	3" (frame base)	3" (motor pivot)	3.25" (frame base)	3" (frame center)	1.5" (anti-tip wheels)
	19.75"	24.5"	35.75"	19.75"	26.5"
foot riggings	35.6"	40"	36"	35.6"	38.25"
	26.25"	24"	25"	26.25"	29.25"
	<b>W:</b> 10-24" <b>D:</b> 10-24"	<b>N/A</b>	<b>W:</b> 12-22" <b>D:</b> 12-22"	<b>W:</b> 20-32" <b>D:</b> 20-28"	<b>W:</b> 20-32" <b>D:</b> 20-28"
0-22"	<b>W:</b> 12-24" <b>D:</b> 12-24"	<b>W:</b> 12-22" <b>D:</b> 12-22"	<b>W:</b> 12-22" <b>D:</b> 12-22"	N/A	N/A
12-22"	<b>W:</b> 12-24" <b>D:</b> 14-24"	<b>W:</b> 12-22" <b>D:</b> 12-22"	<b>W:</b> 12-22" <sup>3</sup> <b>D:</b> 12-22" <sup>3</sup>	N/A	N/A
12-22"	<b>W:</b> 12-24" <b>D:</b> 12-24"	<b>W:</b> 12-22" <b>D:</b> 12-22"	<b>W:</b> 12-22" <sup>3</sup> <b>D:</b> 12-22" <sup>3</sup>	N/A	N/A
	N/A	N/A	N/A	<b>W:</b> 20-32" <b>D:</b> 20-28"	N/A
	N/A	N/A	N/A	<b>W:</b> 20-32" <b>D:</b> 20-25"	N/A
	N/A	N/A	N/A	<b>W:</b> 20-24" <b>D:</b> 20-25" <sup>10</sup>	N/A
	N/A	N/A	N/A	N/A	<b>W:</b> 20-32" <b>D:</b> 20-28"
	N/A	N/A	N/A	N/A	N/A
	N/A	N/A	N/A	N/A	<b>W:</b> 20-32" <b>D:</b> 20-28"
15-20"	N/A	N/A	<b>W:</b> 15-22" <b>D:</b> 15-20"	N/A	N/A
	N/A	N/A	N/A	N/A	N/A
	16.5"-18.5"	15.25-17.25" <small>NOT AVAILABLE WITH SYNERGY®</small>	16.5-18.1"	17.5-19.5"	16.5-18.5"
	18.125"-20.125"	17.5-19.5"	18-19.6"	18.5-20.5"	18.125-20.125"
	19.375"-21.375"	17.75-19.75"	18.5-20.1"	17.75-19.75"	N/A
	N/A	N/A	16.5"-18.25"	N/A	N/A
	N/A	N/A	N/A	N/A	N/A
	Yes	Yes	Yes	No	No
	Group 24	Group 24: 51.5 lbs. Group 34: 38.5 lbs.	Group 24	Group 24	Group 24
	52 lbs. each	51.5 lbs.	52 lbs.	52 lbs.	52 lbs.
3 EX NE NE+	75A Q-Logic 3 NE+ 100A Q-Logic 3 EX	120A Q-Logic 3 EX	75A Q-Logic 2 NE 75A Q-Logic 2 NE+ 100A Q-Logic 2 EX	120A Q-Logic 3 EX 75A Q-Logic 2 NE 75A Q-Logic 2 NE+	90A VR2 6-key 100A Q-Logic 2 EX
	8A Off-board	8A Off-board	8A Off-board	8A Off-board	8A Off-board
Wheel 6®	4-Pole Mid-Wheel 6®	2-pole, front-wheel drive	4-Pole 6 mph <sup>1</sup> or 4-Pole 8 mph <sup>2</sup>	4-Pole 4.7 mph <sup>1</sup>	4-Pole 4.5 mph
	300 lbs.	300 lbs.	300 lbs.	450 lbs. <sup>15</sup>	600 lbs.
	143.5 lbs.	140 lbs.	129 lbs.	143.5 lbs.	195 lbs.

Weights are measured with the unit unloaded.  
 \* Suspension system, seat-to-floor heights  
 0.25" lower when the unit is loaded.  
 \* Height is not available with power tilt  
 \* Seat-to-floor height is not available with  
 power recline.



**QUANTUM®**



888-570-1113

Beamsville, Ontario

[quantumrehab.com](http://quantumrehab.com)

The information contained herein is correct at the time of publication.  
We reserve the right to alter specifications without prior notice.

CA080918