

Jazzy® EVO 613 <sup>1</sup>	Group 2: Standard (K0822/K0823)
Jazzy® EVO 613 Li <sup>1</sup>	Group 2: Standard (K0822/K0823)

Weight Capacity	300 lbs.
Turning Radius <sup>3</sup>	22"
Overall Width <sup>3</sup>	22"
Overall Length <sup>3</sup>	41" with footrest (35" without)
Maximum Speed <sup>2</sup>	Up to 5.1 mph
Ground Clearance <sup>3</sup>	3"
Front Casters <sup>3</sup>	5"
Drive Wheels <sup>3</sup>	13" solid
Rear Casters <sup>3</sup>	5"
Drivetrain	Two motor, in-line, mid-wheel drive
Braking System	Regenerative and electromechanical
Suspension Type	Active-Trac® (ATX)
High-back Seat Sizes	16"x16-18" 18"x18-20" 20"x18-20"
Seat-to-Floor Range	19.75" - 21.25"
Per Charge Range <sup>3,5</sup>	<b>EVO 613:</b> Up to 14 miles at 300 lbs.; Up to 17.1 miles at 200 lbs. <b>EVO 613 Li:</b> Up to 17.6 miles at 300 lbs.; Up to 21 miles at 200 lbs.
Standard Electronics	40A, PG GC3
Component Weights <sup>6</sup>	Base: 94 lbs. Seat: 48 lbs. (18" seat) U1 Batteries: (2) 23.5 lbs. each; Lithium battery: (1) 18.5 lbs.
Total Weight <sup>7</sup>	<b>EVO 613:</b> 189 lbs. <b>EVO 613 Li:</b> 160.5 lbs.
Battery Requirements	<b>EVO 613:</b> (2) 12 volt, U-1 <b>EVO 613 Li:</b> (1) Lithium iron phosphate
Battery Charger	3.5 amp, off-board
Warranty	Frame: 5-year limited Drivetrain: 13-month limited Electronics: 13-month limited SLA Battery: 13-month limited Lithium Battery: 3-year limited

- Jazzy® Power Chairs are FDA Class II Medical Devices designed to aid individuals with mobility impairments.
- Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tire condition. These specifications can be subject to a variance of (+/- 10%).
- Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or - 3%).
- AGM or gel cell type required.
- Tested in accordance with ANSI/RESNA, WC Vol2, section 4 & ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight capacity.
- Battery weight may vary based on manufacturer.
- Total weight includes base, seat and batteries.

**PRIDE**  
MOBILITY

(US) 800-800-8586 • (Canada) 888-570-1113

[pridemobility.com/Jazzy-Evo-613Li](https://pridemobility.com/Jazzy-Evo-613Li) |

HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. DME providers are responsible for determining appropriate billing codes when submitting for insurance reimbursement.

The information contained herein is correct at time of publication; we reserve the right to alter specifications without prior notice.

NOTE: All specifications subject to change without notice. The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3 %. Special order seats or non-standard dimensions will not be fold-flat and depth-adjustable.

010925

**Jazzy**  
**EVO** 613 SERIES

**THE MOST ADVANCED  
JAZZY® EVER!**



FDA Class II  
Medical Device<sup>1</sup>

**ALWAYS INSIST ON LITHIUM!**  
[pridemobility.com/Jazzy-Evo-613-Li](https://pridemobility.com/Jazzy-Evo-613-Li)

# THE MOST ADVANCED JAZZY® EVER!



## Lithium Iron Phosphate Battery\* (Available on Jazzy® EVO 613 Li)

- Faster charging time
- Sustained power performance
- Lighter weight
- More useable energy
- Longer life

Height, width, depth  
and angle adjustable  
armrests

GC3 controller with  
quick disconnect

Depth-adjustable,  
limited recline,  
memory foam  
high-back seat  
with swivel

Narrow base  
width of 22"

13" drive wheels

Removable  
color-through  
shrouds

Convenient access  
to the controller  
and batteries while  
user remains seated


Mid-Wheel 6® Technology  
provides maximum stability

## Always insist on lithium!

The all new Medicare-reimbursable Jazzy® EVO 613 series is here! The Jazzy® EVO 613 series provides the freedom to navigate small spaces and tight corners indoors with a narrow base width of only 22". The Jazzy® EVO 613 Li uses a lithium iron phosphate battery that provides longer life and more convenient charging.

- Maximum speed up to 5.1 mph
- 13" mid-wheel drive tires
- Narrow base width of only 22"
- Removable matte-finish color-through shrouds
- Access to the controller and batteries without the user needing to exit the chair
- Depth-adjustable, limited recline, memory foam single seat post allows the user to swivel
- Pride's patented Active-Trac® (ATX) Suspension
- Mid-Wheel 6® Technology

## Available Colors

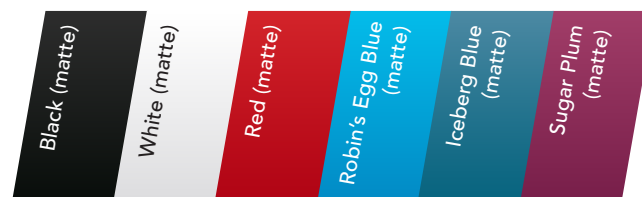

## Accessories

- Cane/crutch holder
- Cell phone holder
- Cup holder
- Flag kit
- Oxygen tank holder
- Rear basket
- Saddlebag
- Walker holder
- Weather cover
- 2.1A XLR USB charger
- 50", 70", 80" and 90" lap belts
- Stealth Products® cushions: Zen, Essence and Glacial

## Options

- Elevated leg rests
- Swing-away leg rests
- Unoccupied transit
- Heel loops
- Residual limb support
- Depth-adjustable solid seat pan, limited recline, high-back with swivel

## Military Patches



\* The Jazzy EVO 613 Li (with lithium batteries) is not airline approved.

FDA Class II Medical Device<sup>1</sup>