

**Jazzy®**  
1170XL PLUS

## COMBINING TIGHT-QUARTER MANEUVERABILITY WITH EXCEPTIONAL OUTDOOR CAPABILITY

### Jazzy® 1170XL Plus

The Jazzy® 1170XL Plus is a stylish, high performance power chair that combines excellent indoor maneuverability with a top speed up to 5 mph and a weight capacity of up to 500 lbs. It is equipped standard with Active-Trac® Suspension which provides for great outdoor performance and rough terrain handling. The Jazzy 1170XL Plus is the ultimate power chair for the active outdoor user who still requires tight-quarter indoor maneuverability.



**Jazzy® 1170XL Plus**  
shown in Viper Blue with gray  
contoured high-back seat, foot  
platform and PG VSI controller.



The Jazzy® 1170XL series models include:

- Jazzy 1170XL Plus 2VHD-C
- Jazzy 1170XL Plus 2VHD-SS
- Jazzy 1170XL Plus 3VHD-C
- Jazzy 1170XL Plus 3VHD-SS

## FEATURES

- 500 lbs. weight capacity
- Powerful, patented, two-motor mid-wheel drive design provides superb power and acceleration
- Active-Trac® and rear caster suspension for outstanding outdoor performance
- Adjustable front anti-tips to meet your preferences and terrain challenges
- On-board battery charger

## OPTIONS


- Contoured high-back with solid seat pan
- Cup holder
- Weathercover
- Rear basket
- Oxygen tank holder
- Cane/crutch holder
- Walker holder
- Swing-away joystick mount
- Safety flag

The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3 %.

# Jazzy® 1170XL Plus

Weight capacity:	500 lbs.
Maximum speed:	Up to 4.5 mph
Ground clearance:	4"
Turning radius:	23"
Overall length:	46"
Overall width:	30.5"
Wheels:	Drive: 16"
Front anti-tips:	9" solid
Rear casters:	8" solid
Suspension:	Active-Trac®
Drive train:	Two-motor, mid-wheel drive
Braking system:	Intelligent braking (electronic, regenerative disc brakes)
Std. electronics:	70 amp, PG VSI controller
Specialty controls:	No
Battery charger:	On-board, 5 amp (standard)
Per-charge range:	Up to 15 miles
Battery requirements:	(2) 12 volt, deep cycle
Battery size:	Group 24
Transport weight	
Battery weight:	53.5 lbs. each
Base weight:	149.5 lbs.
Std. seat weight:	50 lbs.
Max. seat dimensions:	28" x 26"
Warranty:	Lifetime limited warranty on frame; 2-year limited warranty on electronics; 18-months limited warranty on drive motors

Colors:

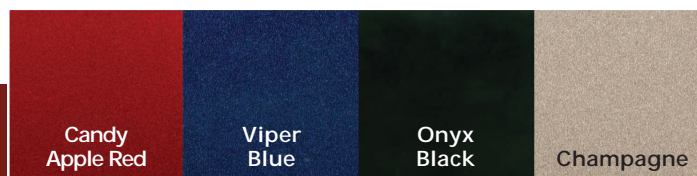

800-800-8586 (US) • 888-570-1113 (Canada) • [www.pridemobility.com](http://www.pridemobility.com)

**Pride®**  
Mobility Products Corp.