

# Jazzy®

## ELITE ES-1




The Jazzy® Elite ES-1 offers a combination of great style, comfort and performance for your active lifestyle.

The in-line, front-wheel drive design provides tight turns around corners for optimal indoor maneuverability in tight spaces. Large front wheels offer extra absorption over obstacles and superb climbing capabilities combine with front anti-tip wheels for transitioning over grade changes easily.

### IN-LINE FRONT-WHEEL DRIVE TECHNOLOGY



PERFORMANCE-MATCHED  
ELECTRONICS

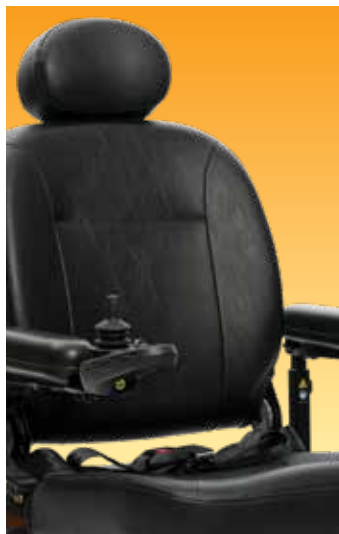

COMFORT HIGH-BACK SEAT  
WITH HEADREST

## AMERICA'S FAVORITE POWER CHAIR

Learn more online.



[gojazzy.com](http://gojazzy.com)

## Specifications

<b>Weight Capacity</b>	300 lbs.
<b>Maximum Speed</b>	Up to 4 mph
<b>Ground Clearance</b>	0.75"
<b>Turning Radius</b>	23.5"
<b>Overall Length</b>	34.75"
<b>Overall Width</b>	22.75"
<b>Drive Wheels</b>	9" solid
<b>Front Wheels</b>	5" solid anti-tips
<b>Rear Wheels</b>	6" solid casters
<b>Suspension</b>	Limited
<b>Drivetrain</b>	Two-motor, in-line, front-wheel drive
<b>Braking System</b>	Regenerative and electro-mechanical
<b>Standard Electronics</b>	40 amp, PG GC 3 controller
<b>Battery Charger</b>	Off-board, 3.5 amp
<b>Per-Charge Range</b>	Up to 12 miles
<b>Battery Requirements</b>	(2) 12 volt, U-1
<b>Battery Weight</b>	25 lbs. each
<b>Base Weight</b>	69 lbs.
<b>Standard Seat Weight</b>	43 lbs.
<b>Standard Seat</b>	16" x 16-18", 18" x 18-20", 20" x 18-20"
<b>Max. Seat Dimensions</b>	20" x 20"
<b>Warranty</b>	5-year limited on frame; 13-month limited on batteries; 13-month limited on electronics; 13-month limited on drive motors

## Colors



Candy Apple Red



Viper Blue

## Models

**Jazzy Elite ES-11S-C**

Code: K0816

\*Additional charges may apply.

## Features and Benefits

- Depth adjustable, limited recline comfort high back seat
- Flip-up foot platform
- Black, flat-free, non-marking tires
- 60" lap belt

## Seating Options\*

- Deluxe contour, limited recline high back available sizes: 16"x16" or 18"x16"

## Options

- Elevating leg rests
- Heel loops
- Residual limb support
- Swing-away leg rests
- Swing-away joystick

## Accessories

- Cane/crutch holder
- Cell phone holder
- Cup holder
- Flag holder
- Lap belt (50", 70", 80", 90")
- Military patches
- Oxygen tank holder
- Rear basket
- Saddle bag
- Walker holder
- Weather cover
- XLR USB charger

## Recommended Lifts

- Backpacker® Series
- Zeus Series
- Outlander Full Platform
- Commander 450



Find a compatible lift for your vehicle online:  
[PRIDEMOBILITY.COM/LIFT-FINDER](http://PRIDEMOBILITY.COM/LIFT-FINDER)

**Pride**  
Mobility Products Corp.

(US) 800-800-8586 • (Canada) 888-570-1113 • [pridemobility.com](http://pridemobility.com)



The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3%. Special order seats or non-standard dimensions will not be fold-flat and depth-adjustable.