

Jazzy[®]

SELECT 6




The Jazzy Select[®] 6 incorporates sporty design and performance features for a truly unique ride experience. Pride's patented Active-Trac[®] Suspension combined with six wheels on the ground and an innovative articulating beam allow for superior performance indoors and out, with added stability over uneven terrain and ramp transitions.

ACTIVE-TRAC[®] SUSPENSION



**JAZZY ARMOR[®] AND
COLOR-THROUGH
SHROUD**

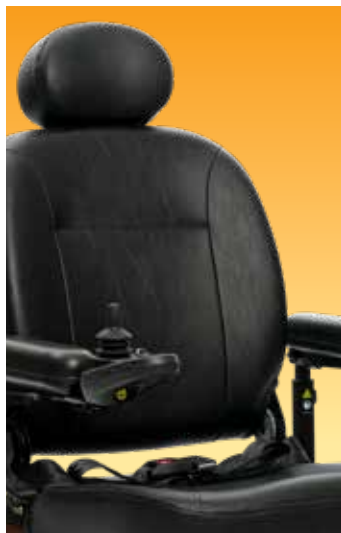

**COMFORT HIGH-BACK SEAT
WITH REPLACEABLE BACK
AND SEAT COVERS**

AMERICA'S FAVORITE POWER CHAIR

Learn more online.



gojazzy.com

Specifications

Weight Capacity	300 lbs.
Maximum Speed	Up to 4 mph
Ground Clearance	1.75"
Turning Radius	21.875"
Overall Length	35.375"
Overall Width	23.5"
Drive Wheels	10"
Front Wheels	6" solid casters
Rear Wheels	6" solid casters
Suspension	Active-Trac® (ATX)
Drivetrain	Two-motor, in-line, mid-wheel drive
Braking System	Regenerative and electro-mechanical
Standard Electronics	40 amp, PG GC 3 controller
Battery Charger	Off-board, 3.5 amp
Per-Charge Range	Up to 15 miles
Battery Requirements	(2) 12 volt, U-1
Battery Weight	25 lbs. each
Base Weight	100 lbs.
Standard Seat Weight	43 lbs. (comfort, high-back)
Standard Seat	16" x 16-18", 18" x 18-20", 20" x 18-20"
Max. Seat Dimensions	20" x 20"
Warranty	5-year limited on frame; 13-month limited on batteries; 13-month limited on electronics; 13-month limited on drive motors

Colors



Jazzy Red



Jazzy Blue

Models

Jazzy Select 6-SS	Code: K0822
Jazzy Select 6-C	Code: K0823

*Additional charges may apply.

Features and Benefits

- Jazzy® Armor® includes shroud and controller guards to protect against daily wear and tear
- Depth adjustable, limited recline comfort high back seat
- Universal, fixed, in-line joystick mount
- Flip-up height adjustable foot platform
- Black, flat-free, non-marking tires
- 60" lap belt

Seating Options*

- Depth adjustable, limited recline comfort high back seat with solid seat pan, available sizes: 16"x16"-18", 18"x18"-20", 20"x18"-20"
- Deluxe contour, limited recline high back available sizes: 18"x16"
- Deluxe contour, limited recline, high back with solid seat pan, available sizes: 16"x16" or 20"x20"
- Power Elevating Seat (note seat does not swivel)
- Synergy® Seat

Options

- Bell shaped foot platform
- Elevating leg rests
- Heel loops
- Residual limb support
- Swing- away footrests
- Swing-away joystick
- VR2 6-key (must and only upgrade with PES option)

Accessories

- Cane/crutch holder
- Cell phone holder
- Cup holder
- Flag holder
- Lap belt (50", 70", 80", 90")
- Military patches
- Oxygen tank holder
- Rear basket
- Saddle bag
- Walker holder
- Weather cover
- XLR USB charger

Recommended Lifts

- Backpacker® Series
- Zeus Series
- Outlander Full Platform
- Commander 450



Find a compatible lift for your vehicle online:
PRIDEMOBILITY.COM/LIFT-FINDER



(US) 800-800-8586 • (Canada) 888-570-1113 • pridemobility.com



The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3%. Special order seats or non-standard dimensions will not be fold-flat and depth-adjustable.