

# Jazzy<sup>®</sup>

## SELECT 6




The Jazzy Select<sup>®</sup> 6 incorporates sporty design and performance features for a truly unique ride experience. Pride's patented Active-Trac<sup>®</sup> Suspension combined with six wheels on the ground and an innovative articulating beam allow for superior performance indoors and out, with added stability over uneven terrain and ramp transitions.

### ACTIVE-TRAC<sup>®</sup> SUSPENSION



**JAZZY ARMOR<sup>®</sup> AND  
COLOR-THROUGH  
SHROUD**

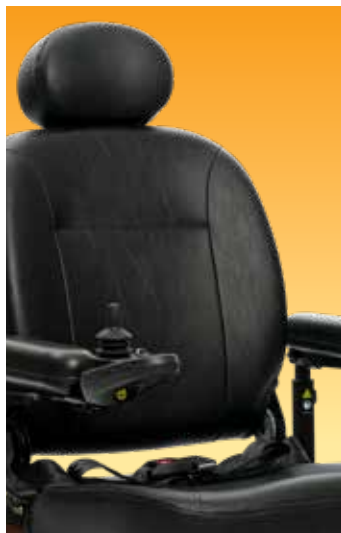

**COMFORT HIGH-BACK SEAT  
WITH REPLACEABLE BACK  
AND SEAT COVERS**

## AMERICA'S FAVORITE POWER CHAIR

Learn more online.



[gojazzy.com](http://gojazzy.com)

# Specifications

<b>Weight Capacity</b>	300 lbs.
<b>Maximum Speed</b>	Up to 4 mph
<b>Ground Clearance</b>	1.75"
<b>Turning Radius</b>	21.875"
<b>Overall Length</b>	35.375"
<b>Overall Width</b>	23.5"
<b>Drive Wheels</b>	10" solid
<b>Front Wheels</b>	6" solid casters
<b>Rear Wheels</b>	6" solid casters
<b>Suspension</b>	Active-Trac® (ATX)
<b>Drivetrain</b>	Two-motor, in-line, mid-wheel drive
<b>Braking System</b>	Regenerative and electro-mechanical
<b>Standard Electronics</b>	40 amp, PG GC 3 controller
<b>Battery Charger</b>	Off-board, 3.5 amp
<b>Per-Charge Range</b>	Up to 15 miles
<b>Battery Requirements</b>	(2) 12 volt, U-1
<b>Battery Weight</b>	25 lbs. each
<b>Base Weight</b>	100 lbs.
<b>Standard Seat Weight</b>	43 lbs. (comfort, high-back)
<b>Standard Seat</b>	16" x 16-18", 18" x 18-20", 20" x 18-20"
<b>Max. Seat Dimensions</b>	20" x 20"
<b>Warranty</b>	5-year limited on frame; 13-month limited on batteries; 13-month limited on electronics; 13-month limited on drive motors

## Colors



Jazzy Red



Jazzy Blue

## Models

**Jazzy Select 6-SS** Code: K0822

**Jazzy Select 6-C** Code: K0823

\*Additional charges may apply.

# Features and Benefits

- Jazzy® Armor® includes shroud and controller guards to protect against daily wear and tear
- Depth adjustable, limited recline comfort high back seat
- Universal, fixed, in-line joystick mount
- Flip-up height adjustable foot platform
- Black, flat-free, non-marking tires
- 60" lap belt

## Seating Options\*

- Depth adjustable, limited recline comfort high back seat with solid seat pan, available sizes: 16"x16"-18", 18"x18"-20", 20"x18"-20"
- Deluxe contour, limited recline high back available sizes: 18"x16"
- Deluxe contour, limited recline, high back with solid seat pan, available sizes: 16"x16" or 20"x20"
- Power Elevating Seat (note seat does not swivel)
- Synergy® Seat

## Options

- Bell shaped foot platform
- Elevating leg rests
- Heel loops
- Residual limb support
- Swing- away footrests
- Swing-away joystick
- VR2 6-key (must and only upgrade with PES option)

## Accessories

- Cane/crutch holder
- Cell phone holder
- Cup holder
- Flag holder
- Lap belt (50", 70", 80", 90")
- Military patches
- Oxygen tank holder
- Rear basket
- Saddle bag
- Walker holder
- Weather cover
- XLR USB charger

## Recommended Lifts

- Backpacker® Series
- Commander Series
- Outlander Series



Find a compatible lift for your vehicle online:  
[PRIDEMOBILITY.COM/LIFT-FINDER](http://PRIDEMOBILITY.COM/LIFT-FINDER)



(US) 800-800-8586 • (Canada) 888-570-1113 • [pridemobility.com](http://pridemobility.com)



The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3%. Special order seats or non-standard dimensions will not be fold-flat and depth-adjustable.