

i L E V E L R E A D Y



**Q6 EDGE® HD**  
shown in *Yellow Blaze*  
with iLevel® technology

# Q6 edge<sup>®</sup> HD

With features like standard 4-pole motors, Mid-Wheel 6<sup>®</sup> Drive Design and ATX Suspension, the Q6 Edge<sup>®</sup> HD is engineered to meet the performance needs of the most active user.

The Q6 Edge HD accepts a complete range of seating and electronics options to deliver exceptional rehab capability at an excellent value, and is available with our iLevel<sup>®</sup> technology.

QUANTUM<sup>®</sup>

# Specifications

<b>DRIVE WHEELS</b>	Standard: 14" Run flat Optional: 14" Pneumatic	
<b>CASTER WHEELS</b>	Front: 6" Solid Rear: 6" Solid	
<b>ANTI-TIP WHEELS</b>	N/A	
<b>MAXIMUM SPEED <sup>2</sup></b>	Up to 4.7 mph	
<b>GROUND CLEARANCE</b>	3" (frame center)	
<b>TURNING RADIUS <sup>3</sup></b>	19.75"	
<b>OVERALL LENGTH</b>	35.6" without foot riggings	
<b>BASE WIDTH</b>	26.25"	
<b>SEATING SIZES <sup>16</sup></b>	<b>TB2 HD</b>	<b>TB3 HD</b>
<b>Synergy<sup>®</sup>/Static Seat</b>	Width: 20-28" Depth: 20-28"	Width: 20-32" Depth: 20-28"
<b>Power Tilt</b>	Width: 20-28" Depth: 20-28"	Width: 20-32" Depth: 20-25"
<b>Lift &amp; Tilt</b>	N/A	Width: 20-24" Depth: 20-25"
<b>SEAT-TO-FLOOR HEIGHTS</b>	<b>TB2 HD</b>	<b>TB3 HD</b>
<b>Synergy<sup>®</sup>/Static Seat <sup>6</sup></b>	16.5-18.5"	17.5-19.5"
<b>Power Tilt <sup>6</sup></b>	18.2-20.9"	18.5-20.5"
<b>Lift &amp; Tilt <sup>6</sup></b>	N/A	17.75-19.75"
<b>MANUAL TILT/RECLINE <sup>7</sup></b>	No	
<b>BATTERY SIZE <sup>4</sup></b>	Group 24	
<b>BATTERY WEIGHT <sup>9</sup></b>	52 lbs.	
<b>AVAILABLE ELECTRONICS <sup>8</sup></b>	75A Q-Logic 2 NE 75A Q-Logic 2 NE+ 100A Q-Logic 2 EX	
<b>BATTERY CHARGER</b>	8A Off-board	
<b>MOTOR PACKAGES</b>	4-Pole 4.7 mph <sup>1</sup>	
<b>WEIGHT CAPACITY <sup>5</sup></b>	450 lbs. <sup>10</sup>	
<b>BASE WEIGHT</b>	143.5 lbs.	
<b>MEDICARE CODE COVERAGE</b>	Group 3: HD (K0850/K0851), Single Power HD (K0858), Multiple Power HD (K0862)	


- 1) Standard
- 2) Speed varies with user weight, terrain type, battery charge, battery condition and tire pressure.
- 3) Without front riggings
- 4) Battery capacity may change when adding power positioning systems.
- 5) Weight capacity is dictated by the selection of seating.
- 6) All measurements +/- 1/8" with 18" x 18" seat, fixed towers, 300 lbs. weight capacity and solid tires (add 3/8" with pneumatic tires). Widths may increase seat-to-floor height. Contact Customer Service for details. Listed measurements are in the lowest position; does not include elevating system unless specified.
- 7) The Synergy<sup>®</sup> Manual Tilt weight capacity is 75 lbs. and maximum dimensions are 14" W x 14" D. Maximum dimensions of the Manual Recline vary by power base and seating system. Call for details.
- 8) All Quantum Rehab<sup>®</sup> units can be ordered with specialty controls
- 9) Battery weight may vary +/- 3 lbs. based on manufacturer.
- 10) Weight capacity is reduced to 400 lbs. with the iLevel<sup>®</sup> power adjustable seat height
- 11) A minimum seating size of 20"x20" is required with iLevel<sup>®</sup>. iLevel<sup>®</sup> maximum seating size is 24" x 25".

# Features

- Mid-Wheel 6<sup>®</sup> allows six wheels on the ground for maximum stability
- Compatible with TRU-Balance<sup>®</sup> 3 HD Power Positioning
- Low impact OMNI-Casters (nylon, spherical-shaped casters) on front and rear prevent wheel hang-ups
- ATX Suspension (Active-Trac<sup>®</sup> with extra stability) incorporates front OMNI-Casters and semi-independent rear caster beam for enhanced performance over more varied terrain
- iLevel<sup>®</sup> option enables up to 10" of seat elevation while at walking speed (up to 3.5 mph)<sup>10,11</sup>

# Options

- iLevel<sup>®</sup> power adjustable seat height<sup>10,11</sup>
- TRU-Balance<sup>®</sup> 3 (450 lbs.) HD Tilt, HD Recline, and HD Tilt & Recline
- Q-Logic 2 Drive Control System
- USB mobile device charger
- Swing-away joystick mount
- LED fender lights (standard with iLevel<sup>®</sup>)

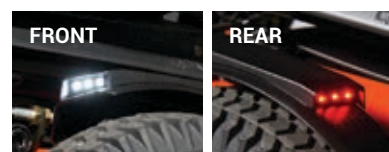

# Colors

SEAT BACKS for units without TRU-Balance<sup>®</sup> 3 are only available in Back in Black.

