



QUANTUM®

Dedicated to revolutionizing the
power chair experience.

quantumrehab.com



T H E 4 F R O N T O F A N E W E R A

4FRONT™

The 4Front™ is a quiet, more responsive front-wheel drive power chair that features automotive-grade suspension with unprecedented comfort and ride quality. The 4Front's CASE (Caster Angle Sensor Encoder) technology and SRS (Smooth Ride Suspension) provide stability for smooth obstacle transitions and unmatched driving performance. The 4Front™ features standard speeds of up to 6 mph in the lowered position, and speeds up to 3.2 mph while elevated at 10" with our optional safe seat elevation. Standard LED fender lights allow you to see and be seen.

QUANTUM® 4FRONT™ shown in *Pretty in Pink* with 10" Safe Seat Elevation and Q-Logic 3 Drive Controls




"The 4Front is by far the most advanced, comfortable and maneuverable chair I have ever driven."

STEWART LUNDY



CONQUER

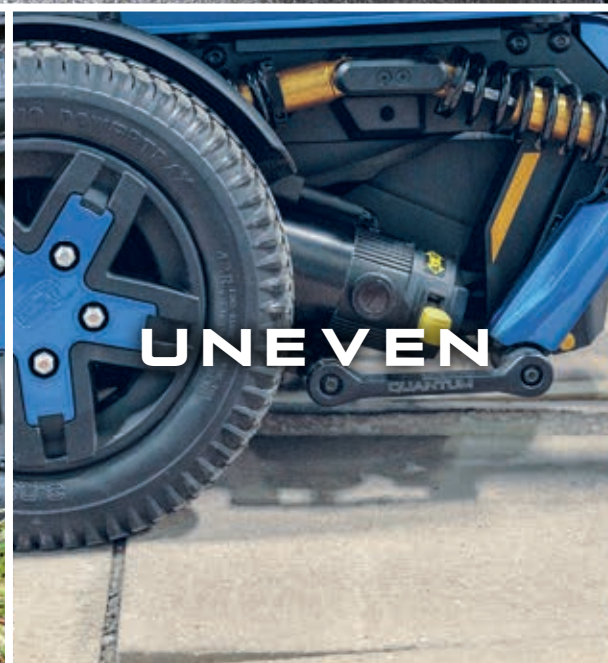

UNEVEN



TERRAIN



iLevel® Power Chairs feature a revolutionary new power chair design for a better experience. iLevel improves daily functionality and quality of life by allowing operation of the power chair with the seat elevated while at walking speed (up to 3.5 mph).

12" LIFT
AT
3.5 MPH

Maximum lift of 10" on Q6 Edge® HD

12" of power adjustable seat height while at walking speed (10" on Q6 Edge® HD)

iLevel's patented Extra Stability Technology™ enhances safety for transfers, reach, and numerous other activities of daily living

Up to 17 exciting color options

Retrofittable on the Q6 Edge® 2.0

300 lbs. weight capacity on Q6 Edge® 2.0 and Q6 Edge® Z

400 lbs. weight capacity on Q6 Edge® HD

Extra Stability Technology™ for a smooth and stable ride

Front and rear LED fender lights are standard on the Q6 Edge® 2.0, Q6 Edge® Z and Q6 Edge® HD with iLevel®





Q6 Edge® 2.0 and Q6 Edge 2.0 X

With features like standard 6 mph motors and advanced responsiveness, the Q6 Edge® 2.0 and 2.0 X are engineered to meet the performance needs of the most active user. The Q6 Edge 2.0 and 2.0 X accept a complete range of seating and electronics options to deliver exceptional rehab capability at an excellent value, and are both available with iLevel® technology.

FEATURES

- User weight capacity of 300 lbs.
- Mid-Wheel 6® allows six wheels on the ground for stability
- Standard 6 mph motors offer advanced performance, precise control and responsiveness via innovative single-stage gearbox
- iLevel® option available with up to 12" of power adjustable seat height at walking speed (3.5 mph)
- Available with Q-Logic 3 Advanced Drive Control System
- Up to 17 exciting color options
- Available with WC-19 crash-test approved occupied transit kit option on tilt & recline and iLevel®
- Available with Accu-Trac Advanced Tracking Technology
- Optional switch-operated front and rear LED fender lights (standard on iLevel®)
- Optional 8" casters

Q6 Edge® Z

The Q6 Edge® Z utilizes powerful 4-pole motors for maximum torque and Mid-Wheel 6® design for aggressive outdoor performance and precision indoor maneuverability. The Q6 Edge® Z accepts a broad range of exclusive seating and electronics options for unmatched adaptability and rehab capability, and is available with our iLevel® technology.

FEATURES

- User weight capacity of 300 lbs.
- Mid-Wheel 6® allows six wheels on the ground for stability
- 8 mph motor package option available
- iLevel® option available with up to 12" of power adjustable seat height at walking speed (3.5 mph)
- 13 exciting color options
- Switch-operated front and rear LED fender lights (standard)
- Available with Accu-Trac Advanced Tracking Technology



REAR-WHEEL DRIVE

power bases

Quantum® Rival®

The Rival® rear-wheel drive power chair features innovative design that provides exceptional outdoor performance and handling while delivering outstanding tight space maneuverability. The Rival accepts a complete range of rehab seating and electronics options.

FEATURES

- 12 exciting color options
- Compact size for outstanding maneuvering in tight spaces
- Maximum speed up to 8 mph with high-speed motor package
- High-performance, high efficiency motor package
- Accepts a wide range of seating and positioning systems
- Two suspension packages
- Compatible with TRU-Balance® 3 Power Positioning Systems
- Sliding battery tray for easy access



BARIATRIC

power bases

Q6 Edge® HD

With features like standard 4-pole motors, Mid-Wheel 6® Drive Design and ATX Suspension, the Q6 Edge® HD is engineered to meet the performance needs of the most active user. The Q6 Edge HD accepts a complete range of seating and electronics options to deliver exceptional rehab capability at an excellent value.

FEATURES

- 450 lbs. weight capacity (400 lbs. with iLevel® option)
- iLevel® option available with up to 10" of power adjustable seat height at walking speed (3.5 mph)
- 13 exciting color options
- Mid-Wheel 6® for maximum stability
- Compatible with TRU-Balance® 3 HD Power Positioning
- ATX Suspension (Active-Trac® with extra stability) incorporates front OMNI-Casters and semi-independent rear caster beam for enhanced performance over more varied terrain
- Optional switch-operated front and rear LED fender lights (standard on iLevel®)



Quantum® 1450

The Quantum® 1450 front-wheel drive bariatric power base features innovative design that provides exceptional outdoor performance while delivering outstanding tight space maneuverability. The 1450 accepts a complete range of rehab seating and electronics options.

FEATURES

- 600 lbs. weight capacity
- Max. speed up to 4.5 mph
- High-performance, high-efficiency motor package
- Compatible with TRU-Balance® Power Positioning and a wide range of other seating and positioning systems



PEDIATRIC

power bases

Quantum® 4Front™

The 4Front™ features standard speeds of up to 6 mph in the lowered position, and speeds up to 3.2 mph while elevated at 10" with our optional safe seat elevation. Standard LED fender lights allow kids to see and be seen.

FEATURES

- CASE (Caster Angle Sensor Encoder) technology
- Unmatched driving performance and ride quality
- Advanced stability for smooth obstacle transitions
- Superior, automotive-grade SRS (Smooth Ride Suspension)
- Available with optional 10" of power adjustable safe seat elevation at speeds of up to 3.2 mph
- Q-Logic 3 EX Advanced Drive Controls
- LED front and rear fender lights
- Quiet ride technology
- Bold aesthetics with 14 color options
- Available with WC-19 crash-test approved, occupied transit option with all configurations



Q6 Edge® 2.0

With features like standard 6 mph motors and advanced responsiveness, the Q6 Edge® 2.0 is engineered to meet the needs of the most active child. The Q6 Edge® 2.0 is available with our iLevel® technology and features bold aesthetics with 13 color options to fit any child's personality. It accepts a complete range of seating and electronics options to deliver exceptional rehab capability at an excellent value.

FEATURES

- iLevel® available for seat sizes as small as 12" x 12" (SOF only)
- Standard 6 mph motors offer advanced performance, precise control and responsiveness via innovative single-stage gearbox
- Bold aesthetics with 17 color options
- iLevel® option available with up to 12" of seat elevation at walking speed (up to 3.5 mph)
- iLevel® is retrofittable on the Q6 Edge® 2.0
- Available with Q-Logic 3 Advanced Drive Control System
- Available with WC-19 crash-test approved occupied transit kit option on tilt & recline and iLevel
- Front and rear LED fender lights (standard)
- Available with Accu-Trac Advanced Tracking
- Optional 8" casters



TRU-BALANCE® 3

power positioning systems



Power Recline

- 85° - 175° of recline
- 300 lbs. weight capacity
- Does not increase seat-to-floor height
- Standard shear reduction and cantilever armrests
- Available in up to 13 back shroud colors
- Available with TRU-Balance® 3 Power Tilt and TRU-Balance® 3 Power Adjustable Seat Lift



Power Tilt

- 50° of tilt
- Lower standard seat-to-floor height
- Unmatched adjustment and durability
- Maintains center of gravity shift
- Available in up to 13 back shroud colors



Power Adjustable Seat Lift with iLevel®

- iLevel® technology enables the unit to drive up to 3.5 mph while fully elevated
- Offers 12" of power adjustable seat lift (10" maximum lift on Q6 Edge® HD)
- Available in up to 13 back shroud colors
- Available on the Q6 Edge® 2.0, Q6 Edge® Z, and Q6 Edge® HD



Power Articulating Foot Platform

- Broad articulation range for enhanced adjustment
- 450 lbs. weight capacity
- Unmatched adjustment, durability and strength
- Multiple calf pads, calf panel and footplate options



Heavy Duty Power Tilt & Recline

- 450 lbs. weight capacity
- 50° of tilt
- 85° - 175° of recline
- Available in up to 13 back shroud colors



Complimentary Cup Holder



TRU-Balance® 3 Power Positioning Systems

- Make fine adjustments while client is in the chair
- Order your seat width independent of your back width at no additional charge
- Internally or externally rotate cantilever armrests for optimal positioning
- Easily retrofit power tilt and power recline
- iLevel® is retrofittable on the Q6 Edge® 2.0
- Easy back height and depth adjustment
- Easy center of gravity adjustment
- Complimentary cup holder
- Premium accessories available
- Optional attendant control
- Tool kit included for quick adjustments



Please Note: Some options may not be available on all power bases. Ask Quantum Customer Service for details.



Q6 Edge® 2.0 shown
with Synergy® Seating
System

Synergy® Seating System

The Quantum® Synergy® Seating System offers a clean, streamlined appearance and unparalleled ease of adjustment. The seating system's unique design allows for height, width, depth and back angle adjustments to be made quickly and easily in the field. The Synergy Seat interfaces with most types of seat and back cushions, and is designed to be the most user-friendly system available today.

FEATURES

- Easily adjustable and versatile design that can grow with the client
- Easy customization
- Several width/depth adjustment packages ranging from 10" to 24"
- Accepts a wide variety of after-market seats and backs
- Heavy Duty option provides a weight capacity of 600 lbs. and seating up to 32" x 28" (model specific)

ADI Q-Back

The ADI Q-Back from Stealth Products, Inc. is available on TRU-Balance® 3 Seating Systems and is offered in three configurations:

- Active Contour ADI Q-Back
- Deep Contour ADI Q-Back
- Extra Deep Contour ADI Q-Back

Available sizes for tilt and static seating:

Heights: 19"-22"
Widths: 16", 18", 20"


Available sizes for power recline seating:

Heights: 22"-25"
Widths: 16", 18", 20"

Note: Additional back sizes are available with Special Order.

TRU-Comfort 2 Back and Seat Cushion

- Adjustable positioning and pressure management
- Compatible with both TRU-Balance® 3 Power Positioning Systems and the Static Seat Frame System
- Available with four-way stretch, recovery, fluid-resistant, breathable cover
- Additional 1" of viscoelastic foam for maximum comfort and pressure management

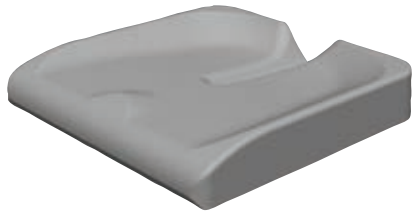


Synergy® Simplicity

Tremendous comfort
and moderate pressure
management

FEATURES

- High-density polyurethane foam construction
- Exclusive wrap-and-lift design for pressure management
- Four-way stretch, fluid resistant upper cover
- Coccyx and sacral relief

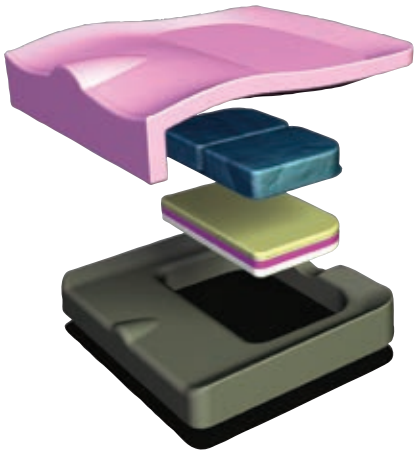


Synergy® Spectrum®

Customized comfort and superb
positioning and pressure
management

FEATURES

- Extra deep well
- Choice of twin-cell lightweight fluid or viscoelastic foam insert
- Viscoelastic top layer
- 2-cover system to prevent shearing with interior liner
- 650 lbs. weight capacity with viscoelastic foam insert

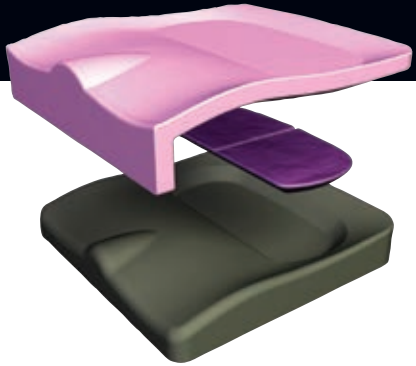


Synergy® Solution®

Three-component design for
pressure management, positioning
and long-lasting comfort

FEATURES

- Exclusive wrap-and-lift design for pressure management
- 2-cover system to prevent shearing with interior liner
- Twin-cell gel insert
- Viscoelastic top layer

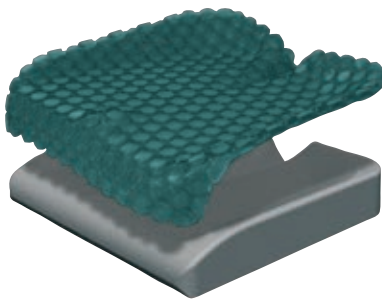

Synergy® Solution® 1

Innovative, two-layered design for
exceptional pressure management
and positioning

FEATURES

- High-density polyurethane foam construction
- Coccyx and sacral relief
- Soft gel top layer for skin protection
- Four-way stretch, fluid resistant upper cover



Active Contour
ADI Q-Back



Deep Contour
ADI Q-Back



Extra Deep Contour
ADI Q-Back



TRU-Comfort 2
Seat Cushion



TRU-Comfort 2
Back Cushion

Featuring
Bluetooth®
(standard)




Q-Logic 3 EX
Joystick Controller



Q-Logic 3 EX
Enhanced Display

Q-Logic 3 Drive Control System

The Q-Logic 3 Advanced Drive Control System from Quantum Rehab® is like no other power wheelchair electronics system. Designed from clinician and consumer feedback, Q-Logic 3 is user-friendly and highly customizable. This adaptable, expandable control system can increase and maximize independence.

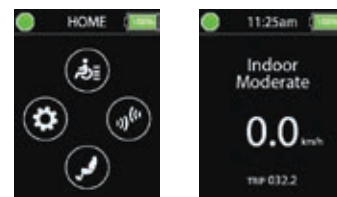

Home Screen

Drive Screen

FEATURES

- Wireless Bluetooth programming through Windows devices
- Bluetooth access to computer, tablet, and phone functions
- iAccess programmable push button and toggle options for customized seating access
- Simple firmware updates
- Capable of utilizing multiple input devices
- Preset drive settings
- Preset color background options
- Clinic Mode:
 - Quick setup of multiple input devices during client evaluations
 - Automatically assigned drive modes to input devices

Q-Logic 3 EX Enhanced Display

The optional Q-Logic 3 EX Enhanced Display serves as a separate screen used with specialty drive controls to display system data. With the Enhanced Display's built-in environmental controls, Q-Logic 3 can operate a wide range of appliances that use infrared (IR) remote controls, such as televisions, and serves as a Bluetooth® computer mouse emulator.



Optional Modules



Q-Logic 3 SAJ



Q-Logic 3 iAccess



Q-Logic 3 Specialty Controls
Interface Module



Non-Expandable



Q-Logic 2
NE and NE+



PG VR2 (4-key)



PG VR2 (6-key)



PG VR2 (9-key)

Expandable (Q-Logic 3)



Q-Logic 3 EX
Joystick



Q-Logic 3 EX
Enhanced
Display



Q-Logic 3 EX
SAJ



Q-Logic 3
Specialty
Controls
Interface
Module

Expandable (Q-Logic 2)



Q-Logic 2 EX
Joystick



Q-Logic 2 EX
Enhanced
Display




Q-Logic 2
EX Stand Alone
Joystick /
Attendant
Joystick



Q-Logic 2
Specialty
Controls
Interface
Module

Multi-Switch Systems

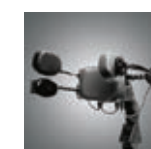

Stealth® i-Drive®
Tri-Array Series,
3-Switch



Stealth® Ultra
3-Switch Head
Array


Proportional Drive Controls



Mo-Vis™
All-Round
Joystick



Mo-Vis™
All-Round Light
Joystick

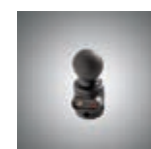

Mo-Vis™
Multi Joystick



Mo-Vis™
Micro Joystick



Stealth®
Precision Mini
Proportional
Joystick



Stealth®
Mushroom
Joystick

Single Switch



Stealth®
Egg Switch



Mo-Vis™
Twister Basic
Button



Mo-Vis™
Twister Satellite
Button

Sip-N-Puff



Stealth®
Sip-N-Puff


Misc. Options



Joystick
Handles


		Edge 2.0	Edge 2.0 X	4Front	Edge Z	Edge HD	Rival™	J6	Q1450
	Sapphire Blue								●
	Deep Purple	●	●	●	●	●	●		
	Midnight Blue	●	●	●	●	●	●		
	Viper Blue							●	●
	Electric Blue	●	●	●	●	●	●		
	True Blue								●
	Emerald Green								●
★	Military Green	●		●	●	●			
★	Night Train	●		●					
	Onyx Black								●
	Back in Black	●	●	●	●	●	●		
	Grey Street	●	●	●	●	●	●		
	Silver								●
	Silver Bullet	●	●	●	●	●	●		
	White	●	●	●	●	●	●		
	Lime Light	●	●	●	●	●	●		
	Yellow Blaze	●	●	●	●	●	●		
	Orion Yellow								●
	Sunburst Orange								●
	Orange Crush	●	●	●	●	●	●		
	Candy Apple Red							●	●
	Red Corvette	●	●	●	●	●	●		
	Pearl Pink								●
	Pretty in Pink	●	●	●	●	●	●		
★	Carbon Fiber	●							
★	Digital Camo	●							
★	Faded Glory	●							

★ New Colors

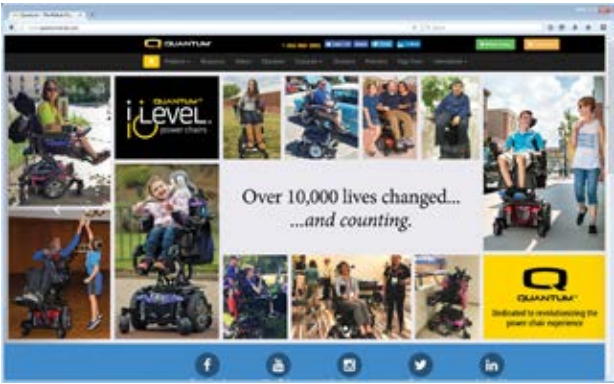


quantumrehab.com

The ultimate resource for everything you need to know about Quantum® products.

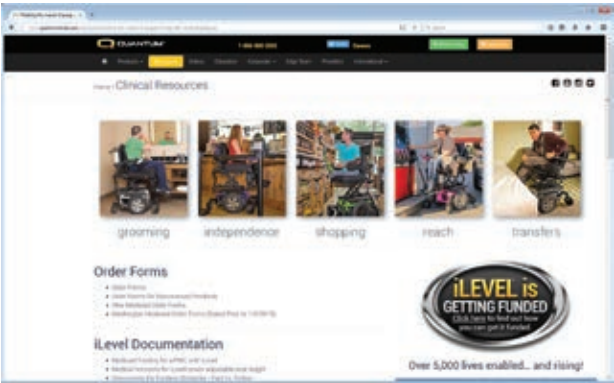


Clinical Resources

A special section of quantumrehab.com for clinicians, featuring order forms, clinical references, and reimbursement information.

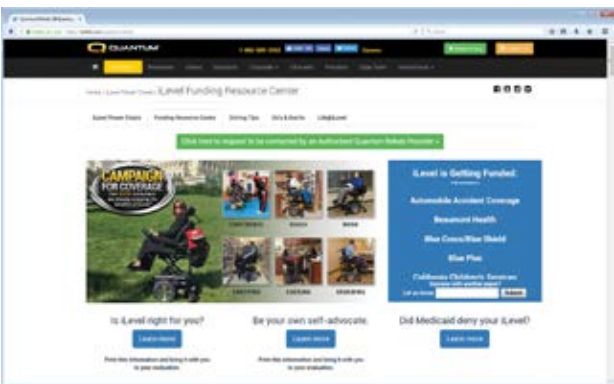

Follow us...

Stay up to date on the latest Quantum® events and product releases through social media.

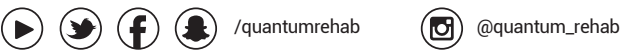

stealthproducts.com

Stealth Products® launched in 1999, in response to a need for versatile, high-quality head and neck positioning systems, as well as switch site location hardware. Their products include: head and body support/positioning, joystick mounts, trays, seating, and pediatric mobility products.

Quantum Rehab® Professional App

The Quantum Rehab® Professional App is available for Apple and Android smart phones and tablets. Look for the app in the Apple App Store or Google Play links using your smart device to download the app directly.

Sales Support (866-800-2002)

Assistance with design, ordering, pricing and quality assurance.

Custom Group (877-465-2999)

Providing innovative custom solutions and specialized product modifications – your “off the order form” solution team.

Technical Support (ext. 5200)

Assistance with troubleshooting and parts order placement.

Product Planning and Reimbursement Center

Assistance with understanding Medicare codes and coverage criteria, product planning, provision guidelines, claims coding, pre-screens and assistance in regard to high-end rehab clients.

Set to Spec Service

Set to Spec Service from Quantum Rehab® is designed to assist you in delivering Q6 Power Chairs to your clients. We will pre-set, to your specifications, seat-to-floor heights, seating center of gravity, armrest heights, front riggings, thoracic and/or thigh laterals.



	Q6 EDGE® 2.0 <small>AVAILABLE WITH iLEVEL®</small>	4FRONT™ <small>AVAILABLE WITH 10" LIFT</small>	Q6 EDGE®	Q6 EDGE® 2.0X	Q6 EDGE® Z <small>AVAILABLE WITH iLEVEL®</small>	QUANTUM® 610	Q6 EDGE® HD <small>AVAILABLE WITH iLEVEL® (400 LBS. WT. CAPACITY)</small>	QUANTUM® RIVAL®	J6	QUANTUM® 1450
DRIVE CONFIGURATION	Mid-wheel Drive	Front-wheel Drive	Mid-wheel Drive	Mid-wheel Drive	Mid-wheel Drive	Mid-wheel Drive	Mid-wheel Drive	Rear-wheel Drive	Mid-wheel Drive	Front-wheel Drive
DRIVE WHEELS	14" Solid Optional: 14" Pneumatic; 14" Knobby	14" Flat Free	14" Solid Optional: 14" Pneumatic	14" Solid Optional: 14" Pneumatic; 14" Knobby	14" Run Flat Optional: 14" Pneumatic	10" Solid Optional: 10" Pneumatic	14" Run Flat Optional: 14" Pneumatic	14" Solid Optional: 14" Pneumatic	10" Solid Optional: 10" Pneumatic	14" Solid
CASTER WHEELS	Front: 6" Solid (Optional: 8") Rear: 6" Solid (Optional: 8")	Rear: 9" Flat Free	Front: 6" Solid Rear: 6" Solid	Front: 6" Solid (Optional: 8") Rear: 6" Solid (Optional: 8")	Front: 6" Solid Rear: 6" Solid	Front: 5" Solid Rear: 6" Solid	Front: 6" Solid Rear: 6" Solid	Front: 9" Pneumatic	Front: 5" Solid Rear: 6" Solid	9" Solid
ANTI-TIP WHEELS	N/A	Front: 4"	N/A	N/A	N/A	N/A	N/A	N/A	N/A	4" Solid
MAXIMUM SPEED ⁴	Up to 6 mph	Up to 6 mph	Up to 5 mph (Optional: 6 mph)	Up to 6 mph	Up to 6 mph (Optional: 8 mph) ¹³	Up to 4 mph	Up to 4.7 mph	Up to 6 mph (Optional: 8 mph)	Up to 4 mph	Up to 4.5 mph
GROUND CLEARANCE	2.7" (battery tray)	3" (motor pivot)	3" (frame center)	2.7" (battery tray)	3" (frame base)	2.5" (frame)	3" (frame center)	3.25" (frame base)	2.5" (frame)	1.5" (anti-tip wheels)
TURNING RADIUS ⁵	6" casters: 20.5" 8" casters: 23.25"	24.5"	20"	6" casters: 20.5" 8" casters: 23.25"	19.75"	22"	19.75"	35.75"	22"	26.5"
OVERALL LENGTH ¹¹	6" casters: 35.5" 8" casters: 40.48"	40"	35.5"	6" casters: 35.5" 8" casters: 40.48"	35.6"	34.74"	35.6"	36"	34.74"	38.25"
BASE WIDTH	24.4"	24"	24"	24.4"	26.25"	23.19"	26.25"	25"	23.19"	29.25"
Seating Sizes:	TRU-Balance® 3	TRU-Balance® 3	TRU-Balance® 3	TRU-Balance® 3	TRU-Balance® 3	TRU-Balance® 2	TRU-Balance® 2 HD TRU-Balance® 3 HD	TRU-Balance® 3	TRU-Balance® 2	TRU-Balance® 2 HD Tilt
SYNERGY®	W: 10-22" D: 10-22"	N/A	W: 10-22" D: 10-22"	W: 10-22" D: 10-22"	W: 10-24" D: 10-24"	W: 10-20" D: 10-20"	W: 20-28" D: 20-28" W: 20-32" D: 20-28"	W: 12-22" D: 12-22"	W: 10-20" D: 10-20"	W: 20-32" D: 20-28"
STATIC	W: 10-22" D: 10-22"	W: 12-22" D: 12-22"	W: 10-22" D: 10-22"	W: 10-22" D: 10-22"	W: 12-24" D: 12-24"	W: 10-20" D: 10-20"	W: 20-28" D: 20-28" W: 20-32" D: 20-28"	W: 12-22" D: 12-22"	W: 10-20" D: 10-20"	W: 20-32" D: 20-28"
POWER TILT	W: 12-22" D: 12-22"	W: 12-22" D: 12-22"	W: 12-22" D: 12-22"	W: 12-22" D: 12-22"	W: 12-24" D: 14-24"	W: 14-20" D: 14-20"	W: 20-28" D: 20-28" W: 20-32" D: 20-25"	W: 12-22" ³ D: 12-22" ³	W: 14-20" D: 14-20"	W: 20-32" D: 20-28"
LIFT & TILT	W: 12-22" D: 12-22"	W: 12-22" D: 12-22"	W: 12-22" D: 12-22"	W: 16-22" D: 14-22"	W: 12-24" D: 12-24"	W: 16-20" D: 16-20"	N/A W: 20-24" D: 20-25" ¹⁰	W: 12-22" ³ D: 12-22" ³	N/A	N/A
Seat-to-Floor Heights ¹⁶ :	TRU-Balance® 3	TRU-Balance® 3	TRU-Balance® 3	TRU-Balance® 3	TRU-Balance® 3	TRU-Balance® 2	TRU-Balance® 2 HD TRU-Balance® 3 HD	TRU-Balance® 3	TRU-Balance® 2	TRU-Balance® 2 HD Tilt
SYNERGY® / STATIC ⁸	16.5-18.5"	15.25-17.25" <small>NOT AVAILABLE WITH SYNERGY®</small>	16.5-18.5"	16.5-18.5"	16.5"-18.5"	16.5-17.5"	16.5-18.5" 17.5-19.5"	16.5-18.1"	16.5-17.5"	16.5-18.5"
POWER TILT ⁸	17.5-19.5"	17.5-19.5"	17.5-19.5"	17.5-19.5"	18.125"-20.125"	16.875"	18.2-20.9" 18.5-20.5"	18-19.6"	16.875"	18.125-20.125"
LIFT & TILT ⁸	17.5-19.5"	17.75-19.75"	18-20"	17.5"-19.5"	19.375"-21.375"	17.75"	N/A 17.75-19.75"	18.5-20.1"	N/A	N/A
MANUAL RECLINE	Yes	Yes	Yes	Yes	Yes	Yes	No	Yes	Yes	No
BATTERY SIZE ⁶	NF-22	Group 24: 51.5 lbs. Group 34: 38.5 lbs.	NF-22	NF-22	Group 24	U-1	Group 24	Group 24	U-1	Group 24
BATTERY WEIGHT ¹⁴	38 lbs.	51.5 lbs.	38 lbs.	38 lbs.	52 lbs. each	23.4 lbs.	52 lbs.	52 lbs.	23.4 lbs.	52 lbs.
AVAILABLE ELECTRONICS ¹²	120A Q-Logic 3 EX 75A Q-Logic 2 NE 75A Q-Logic 2 NE+	120A Q-Logic 3 EX	75A Q-Logic 2 NE 75A Q-Logic 2 NE+ 100A Q-Logic 2 EX	120A Q-Logic 3 EX	75A Q-Logic 2 NE+ 100A Q-Logic 2 EX	70A PG VR2 70A Q-Logic 2 EX	120A Q-Logic 3 EX 75A Q-Logic 2 NE 75A Q-Logic 2 NE+	75A Q-Logic 2 NE 75A Q-Logic 2 NE+ 100A Q-Logic 2 EX	70A Q-Logic 2 NE 70A Q-Logic 2 NE+ 70A Q-Logic 2 EX	90A VR2 6-key 100A Q-Logic 2 EX
BATTERY CHARGER	8A Off-board	8A Off-board	8A Off-board	8A Off-board	8A Off-board	5A Off-board	8A Off-board	8A Off-board	5A Off-board	8A Off-board
MOTOR PACKAGES	2-Pole Mid-Wheel 6®	2-pole, front-wheel drive	4-Pole 5 mph ¹ or 4-Pole 6 mph ²	2-Pole Mid-Wheel 6®	4-Pole Mid-Wheel 6®	2-Pole Glide 4 mph ¹	4-Pole 4.7 mph ¹	4-Pole 6 mph ¹ or 4-Pole 8 mph ²	2-Pole 4 mph	4-Pole 4.5 mph
WEIGHT CAPACITY ⁷	300 lbs.	300 lbs.	300 lbs.	300 lbs.	300 lbs.	300 lbs. ⁹	450 lbs. ¹⁵	300 lbs.	300 lbs. ⁹	600 lbs.
BASE WEIGHT	138.6 lbs. (6" casters) 148 lbs. (8" casters)	140 lbs.	138 lbs.	138.6 lbs. (6" casters) 148 lbs. (8" casters)	143.5 lbs.	100.2 lbs.	143.5 lbs.	129 lbs.	100.2 lbs.	195 lbs.
MEDICARE CODE COVERAGE - For more information on Medicare code coverage, please consult our Product Planning and Reimbursement Center.										
	Group 3: Standard (K0848/K0849), Single Power (K0856), Multiple Power (K0861)	Group 3: Standard (K0848/K0849), Single Power (K0856), Multiple Power (K0861)	Group 3: Standard (K0848/K0849), Single Power (K0856/K0857), Multiple Power (K0861)	Group 3: Standard (K0848/K0849), Single Power (K0856), Multiple Power (K0861)	Group 3: Standard (K0848/K0849), Single Power (K0856), Multiple Power (K0861)	Group 2: Standard (K0822/K0823), Single Power (K0835/K0836)	Group 3: Heavy Duty (K0850/K0851), Single Power HD (K0858), Multiple Power HD (K0862)	Group 3: Standard (K0848/K0849), Single Power (K0856), Multiple Power (K0861)	Group 2: Standard (K0822/K0823), Single Power (K0835)	Group 3: Very Heavy Duty (K0852/K0853), Single Power Very HD (K0860), Multiple Power Very HD (K0863)

1) Standard

2) Optional

3) Requires HD Suspension

4) Speed varies with user weight, terrain type, battery charge, battery condition and tire pressure.

5) Without front riggings

6) Battery capacity may change when adding power positioning systems.

7) Weight capacity is dictated by the selection of seating.

8) All measurements +/- 1/8" with 18" x 18" seat, fixed towers, 300 lbs. weight capacity and solid tires (add 3/8" with pneumatic tires). Depths may increase seat-to-floor height. Contact Customer Service for details. Listed measurements are in the lowest position; does not include elevating system unless specified.

9) 250 lbs. weight capacity with TRU-Balance® Systems

10) A minimum seating size of 20" x 20" is required with iLevel® on the Q6 Edge® HD; maximum seating size is 24" x 25"

11) Without foot riggings

12) All Quantum units can be ordered with specialty controls

13) iLevel® not available with 8 mph motor package

14) Battery weight may vary +/- 3 lbs. based on manufacturer.

15) Weight capacity is reduced to 400 lbs. with iLevel® power adjustable seat height

16) Seat-to-floor heights are measured with the unit unloaded. Due to the 4Front's suspension system, seat-to-floor heights are approximately 0.25" lower when the unit is loaded.

HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. DME providers are responsible for determining appropriate billing codes when submitting for insurance reimbursement.

The information contained herein is correct at the time of publication. We reserve the right to alter specifications without prior notice.



QUANTUM[®]

(US) 866-800-2002 • (CAN) 888-570-1113

Learn more online at: quantumrehab.com



The information contained herein is correct at the time of publication. We reserve the right to alter specifications without prior notice.

QR-1-4-18