

Colt EXECUTIVE

ADVANCED ENGINEERING FOR POWERFUL PERFORMANCE AND EXCEPTIONAL STYLE

COLT EXECUTIVE

Bristling with high-performance features like full suspension, low-profile tyres and a hydraulic sealed brake system, the **Colt Executive** is a blend of power and precision. The Colt Executive features sleek style and a wealth of standard luxury touches like a wraparound delta tiller, new double A-arm suspension, new terrain gripping pneumatic tyres and LED lighting, making it the ultimate outdoor scooter.



Colt Executive
Shown in Titanium Grey
(Images are for illustration
purposes only)


New Terrain Gripping Tyres

New Double A-Arm Suspension

Colt Executive

Weight capacity:	28 Stone 8 lbs. (181 kg)
Maximum speed ¹ :	Up to 8 mph (12.8 km/h)
Ground clearance ² :	13.33 cm (5.25") at side frame
Turning radius:	210.8 cm (83")
Overall length ² :	142.8 cm (56.25")
Overall width ² :	69.2 cm (27.25")
Front tyres:	10 cm x 33 cm (4" x 13") (pneumatic)
Rear tyres:	10 cm x 33 cm (4" x 13") (pneumatic)
Suspension:	Front and Rear
Total Weight Without Batteries ⁵ :	117 kg (259.5 lbs.)
Heaviest piece when disassembled:	Main frame: 93.21 kg (205.5 lbs.)
Standard seating:	Type: Pillow top captain style Weight: 21.77 kg (48 lbs.) Material: Black Vinyl Dimensions: width - 45.72 cm (18") (overall width) depth - 44.45 cm (17.5") (effective depth) height - 45.72 cm (18")
Seat-to-deck range:	49.5 cm - 53.3 cm (19.5" - 21")
Drive system:	Rear-wheel drive, sealed transaxle, 24-volt DC motor
Dual braking system:	Regenerative and electromechanical
Battery requirements ³ :	Two 12-volt, deep cycle Size: 100 AH (standard) Weight ⁴ : 100 AH: 27.44 kg (60.5 lbs.) each
Battery charger:	8-amp, off-board

Colours:

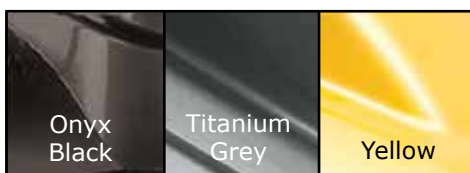

Features

- New Double A-Arm Suspension for optimum handling
- Newly redesigned tiller with LCD console
- New large 33 cm (13") terrain gripping pneumatic tyres for excellent outdoor performance
- Full suspension for optimum handling
- Sealed hydraulic brakes
- Full lighting package including directional signals and lower headlight
- Delta tiller with wraparound handles
- Easily accessible tie down points (for transport of unoccupied scooter)
- Rearview mirror standard
- Front shield standard

Options

- Double crutch holder
- Oxygen tank holder
- Rear basket
- Additional rearview mirror



NEW LCD console with speed sensor (standard)

¹Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tyre condition. This specification can be subject to a variance of +10%,-5%.

²Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or - 3%).

³AGM or gel cell type required.

⁴Battery weight may vary based on manufacturer.

⁵Includes standard seat.

NOTE: All specifications subject to change without notice.

The information contained herein is correct at the time of publication; due to the manufacturer's commitment to constant improvement and development, they reserve the right to alter specifications without prior notice. Theoretical range is calculated under test conditions in accordance with European standards. Actual range in normal use will depend on many factors, including the condition of the vehicle and its batteries, the weight of the driver, correct tyre pressure, ambient temperature, and the gradient and surface of the road or pavement.

32 Wedgwood Road, Bicester, Oxfordshire OX26 4UL
(Tel) 01869 324 600 • (Fax) 01869 323 070 • www.pridemobility.co.uk



DEALER INFO HERE

