

QUANTUM[®]

610

TRUE ALL AROUND PERFORMANCE



Q610 shown with Synergy[®] Seating System,
and PG VR2 controller

The compact Quantum[®] 610 is the ultimate choice for users who demand the superb stability of Mid-Wheel 6[®] Design and tight-quarter manoeuvrability. Highly adaptable, the Quantum 610 accepts a wide range of seating and electronic options while its powerful drivetrain delivers high-performance torque.



FEEL THE POWER

Specifications

DRIVE WHEELS	25.4 cm (10")
OMNI-CASTORS	Front: 12.7 cm (5") solid Rear: 15.24 cm (6") solid
SUSPENSION	ATX Suspension (Active-Trac® with extra stability)
MAXIMUM SPEED ¹	Up to 4 mph (6.44 km/h)
BRAKES	Intelligent braking (electronic, regenerative disc brakes)
GROUND CLEARANCE	6 cm (2.375")
TURNING RADIUS ²	53.98 cm (21.25")
OVERALL SIZE ²	Length: 88.27 cm (34.75") Width: 59 cm (23.25")
SEATING	Synergy® Seat Adjustable Width: 25 cm - 51 cm (10"-20") Adjustable Depth: 25 cm - 51 cm (10"-20") TRU-Comfort Plus Seating ⁵ Adjustable Width: 41 cm-56 cm (16"-20") Adjustable Depth: 41 cm-56 cm (16"-20") TRU-Comfort Plus Back Heights 61 cm (24") or 68.5 cm (27") TRU-Balance® Power Positioning System Adjustable Width: 35.5 cm - 51 cm (14"-20") Adjustable Depth: 35.5 cm - 51 cm (14"-20")
DRIVETRAIN	2-motor, Mid-Wheel 6®
BATTERIES	Two 12 volt, deep cycle, 30 AH
BATTERY CHARGER	5A, off-board (std.)
AVAILABLE ELECTRONICS	70 amp, PG VR2 (standard) 70 amp, Pilot Plus (optional)
WEIGHT CAPACITY	21 Stone 6 lbs. (136 kg)
BASE WEIGHT	51 kg (112 lbs.)
BATTERY WEIGHT ⁴	11 kg (25 lbs.) (each)
WARRANTY	Lifetime limited warranty on frame; 1-year limited warranty on electronics; 18-month warranty on drive control system; 2-year limited warranty on drivetrain


Features

- Mid-Wheel 6® allows six wheels on the ground for maximum stability
- Compatible with TRU-Balance® Power Positioning
- OMNI-Castors (nylon, spherical-shaped castors) on front and rear prevent wheel hang-ups
- Side-mounted, easily accessible freewheel levers
- ATX Suspension (Active-Trac® with extra stability) incorporates front OMNI-Castors for enhanced performance over more varied terrain
- Powerful motors for added power and torque
- 5A, Off-board charger

Options

- Synergy® Seating System
- TRU-Comfort Plus Seating System
- TRU-Balance® Power Positioning System
- Cup holder
- Oxygen bottle holder
- Crutch holder
- Rear basket
- Leg rest options (Please see prescription form)
- Cantilever arm rests
- Synergy® centre mount adjustable foot platform
- Colour shrouds (for TRU-Comfort Plus only)³
- Ottobock arm pads
- Gel arm pads
- Paediatric Tilt
- 12.7 cm (5") post lift (only available with standard version)

Shroud Colors

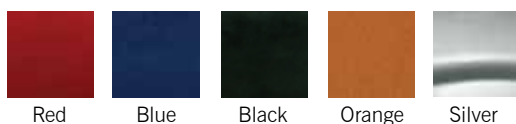

FOOTNOTES:

- 1) Range and speed vary with user weight, terrain type, battery charge, battery condition and tyre pressure.
- 2) Due to manufacturing tolerances and continual product improvement, this specification can be subject to variance of + or - 3 %. Dimensions listed are for power base only. Overall measurements will vary based on seating and accessory selections.
- 3) Black shrouds standard.
- 4) Battery weight may vary based on manufacturer.
- 5) Only available in 16"W x 16"D, 18"W x 18"D or 20"W x 20"D.

The information contained herein is correct at the time of publication; due to the manufacturer's commitment to constant improvement and development, they reserve the right to alter specifications without prior notice. Theoretical range is calculated under test conditions in accordance with European standards. Actual range in normal use will depend on many factors, including the condition of the vehicle and its batteries, the weight of the driver, correct tyre pressure, ambient temperature, and the gradient and surface of the road or pavement.

DEALER INFO HERE



32 Wedgwood Road, Bicester,
Oxfordshire OX26 4UL
(Tel) 01869 324 600
(Fax) 01869 323 070
www.quantumrehab.co.uk

