

# Full Platform Exterior Lifts

Base weight lift	92 lbs. <sup>1</sup>
Mobility device overall weight	380 lbs. <sup>6, 7</sup>
Maximum product width	28"
Maximum travel distance from the ground	20"
Motor	Sealed 12 volt DC
Overall platform length	49.25" (41" usable surface)
Warranty	3-year transferable
Hitch Class <sup>2</sup>	Class II <sup>3</sup> , III <sup>8</sup> , IV <sup>6</sup>
Vehicle platform	Full-size van, Mini-van, Truck, Sedan, SUV and Crossover

<sup>1</sup> Weight does not include any accessories

<sup>2</sup> Vehicle dependent

<sup>3</sup> 190 lbs. maximum product weight with a Class II hitch

<sup>4</sup> Call for details when lifting units with power positioning systems

<sup>5</sup> Optional lockdown arm required. Call for details.

<sup>6</sup> Standard Class IV receiver only

<sup>7</sup> Vehicle requires Class IV hitch with a tongue weight of at least 500 lbs.

<sup>8</sup> 350 lbs. maximum product weight with a Class III hitch

## Can easily lift the following current Pride® Power Chairs<sup>2, 4</sup>

- Jazzy® 600
- Jazzy® 614/614 HD
- Jazzy Select®
- Jazzy Select® 6
- Jazzy Select® 6 Ultra
- Jazzy Select® 14
- Jet™ 3 Ultra

## Can easily lift the following current Pride® Scooters<sup>5</sup>

- Go-Go Elite Traveller®/ LR 3 & 4
- Go-Go Elite Traveller® Plus/ Plus HD 3 & 4
- Go-Go® Ultra X 3 & 4
- Celebrity® X 3 & 4
- Legend® 3 & 4
- Maxima 3 & 4
- Revo® 3 & 4
- Victory® 9 3 & 4
- Victory® 10 3 & 4
- Victory® PS 3 & 4

## Can easily lift the following current Quantum® Power Chairs<sup>2, 4</sup>

- Quantum® 600/600XL
- Quantum® 610
- Quantum® 640
- Quantum® 1103 Ultra
- Quantum® 1121
- Quantum® 6000 XL
- Quantum® 6000Z/6400Z
- Quantum® 6400
- Quantum® R4000/R4400

## Options

- License plate kit
- Optional battery pack with charger
- Swing-away joint
- Extended swing-away joint
- Power chair weather cover
- Scooter weather cover

<sup>1</sup> Weight does not include any accessories

<sup>2</sup> Vehicle dependent

<sup>3</sup> 190 lbs. maximum product weight with a Class II hitch

<sup>4</sup> Call for details when lifting units with power positioning systems


<sup>5</sup> Optional lockdown arm required. Call for details.

<sup>6</sup> Standard Class IV receiver only

<sup>7</sup> Vehicle requires Class IV hitch with a tongue weight of at least 500 lbs.

<sup>8</sup> 350 lbs. maximum product weight with a Class III hitch

# Full Platform Exterior Lifts

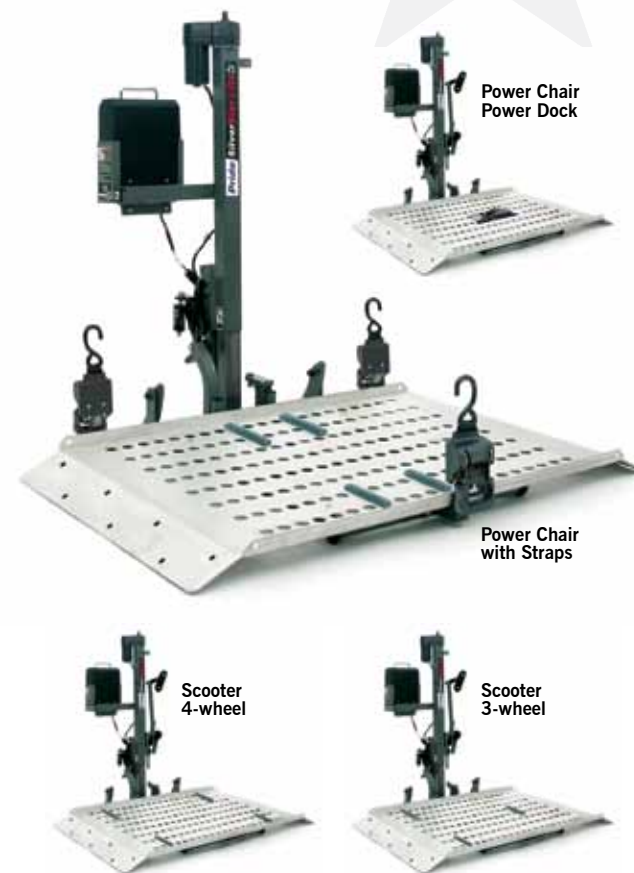

Power Chair  
Power Dock

Power Chair  
with Straps

Scooter  
4-wheel

Scooter  
3-wheel

# Full Platform Exterior Lifts

Silver Star®'s Full Platform exterior lift was designed to provide a durable, easy-to-use system to take your power chair or scooter with you. Simple hitched-based installation requires little to no modification of your vehicle and allows for easy removal.

- 380 lbs. lifting capacity\*
- 3-year transferable warranty
- Drive-on/drive-off platform allows user to enter from either direction
- Easy to use, three-point tie down system with new one-handed, black ratcheting straps
- Lock-down arm available for use with scooters
- Optional Power Dock now available
- Standard adjustable wheel chocks
- Manual override
- Can lift a wide variety of Pride power chairs, scooters and other manufacturers' mobility devices
- Shock-absorbing "cradle effect" protects power chairs from jarring vibrations
- Optional battery pack means no vehicle hardwiring


\*Vehicle requires Class IV hitch with a tongue weight of at least 500 lbs.

Notes: The Full Platform Exterior Lift can lift all of the listed scooters or power chairs with contoured, Synergy® or Specialty seating.

The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice.




Standard adjustable wheel chocks.

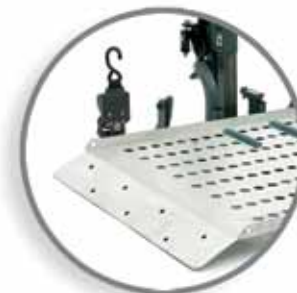

New drive-on/drive-off platform.



New, one-handed, black ratcheting straps.

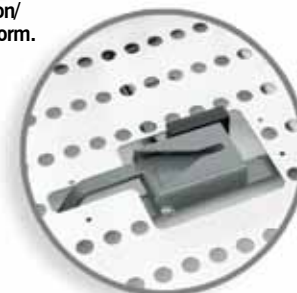

Optional Power Dock now available.

